

(19) World Intellectual Property Organization
International Bureau



PCT



(43) International Publication Date
14 February 2008 (14.02.2008)

(10) International Publication Number
WO 2008/019015 A2

(51) International Patent Classification:
G06F 19/00 (2006.01) **A61M 5/172** (2006.01)

(74) Agent: **BRUESS, Steven, C.**; Merchant & Gould P.c., P.O. Box 2903, Minneapolis, MN 55402-0903 (US).

(21) International Application Number:
PCT/US2007/017123

(81) Designated States (*unless otherwise indicated, for every kind of national protection available*): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BI, BR, BW, BY, BZ, CA, CI, CN, CO, CR, CU, CZ, DE, DK, DM, DO, DZ, EC, EE, EG, ES, FI, GB, GD, GE, GH, GM, GT, HN, HR, HU, ID, IL, IN, IS, JP, KE, KG, KM, KN, KP, KR, KZ, LA, LC, LK, LR, LS, LT, LU, LY, MA, MD, ME, MG, MK, MN, MW, MX, MY, MZ, NA, NG, NI, NO, NZ, OM, PG, PH, PL, PT, RO, RS, RU, SC, SD, SE, SG, SK, SL, SM, SV, SY, TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VC, VN, ZA, ZM, ZW.

(22) International Filing Date: 31 July 2007 (31.07.2007)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:
60/835,927 3 August 2006 (03.08.2006) US
11/702,922 5 February 2007 (05.02.2007) US

(71) Applicant (*for all designated States except US*): **SMITHS MEDICAL MD, INC.** [US/US]; 1265 Grey Fox Road, St. Paul, Minnesota 55112 (US).

(84) Designated States (*unless otherwise indicated, for every kind of regional protection available*): ARIPO (BW, GH, GM, KE, LS, MW, MZ, NA, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HU, IE, IS, IT, LT, LU, LV, MC, MT, NL, PL, PT, RO, SE, SI, SK, TR), OAPI (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).


(72) Inventors; and

(75) Inventors/Applicants (*for US only*): **EVANS, William, James** [GB/US]; 4330 19th Street, San Francisco, CA 94114 (US). **BLOMQUIST, Michael, L.** [US/US]; 11050 Amen Circle, Blaine, MN 55449 (US). **POPE, Rhall, E.** [US/US]; 4825 East Lake Harriet Parkway, Minneapolis, MN 55409 (US).

Published:

— *without international search report and to be republished upon receipt of that report*

(54) Title: INTERFACE FOR MEDICAL INFUSION PUMP

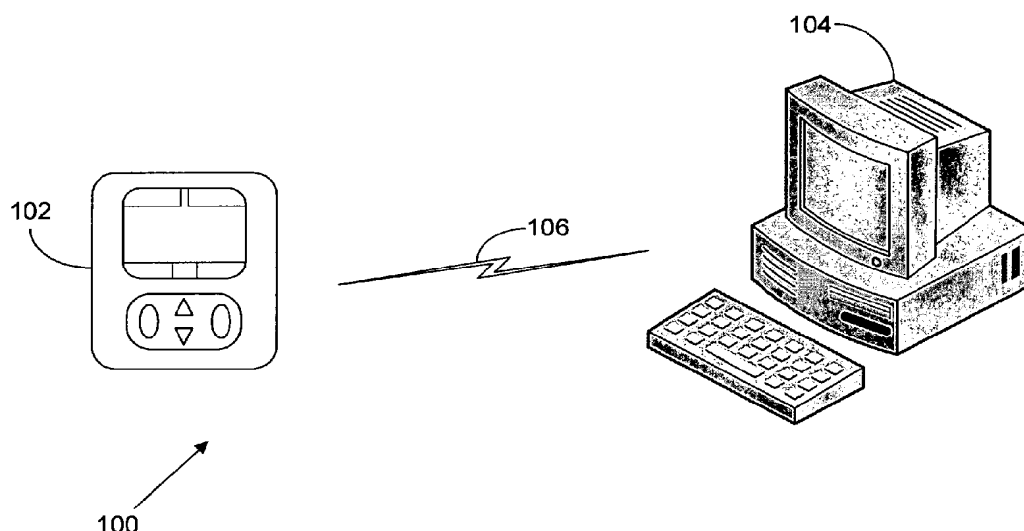

(57) Abstract: An apparatus and method for indicating a current condition of a medical infusion pump is disclosed. The apparatus includes a memory configured to store two or more color codes, each color code corresponding to a condition of the medical infusion pump and a color. The apparatus further includes a monitor. The apparatus also includes a programmable circuit in electrical communication with the memory and the monitor, the programmable circuit programmed to generate a graphical user interface displaying a color upon detection of a condition. The method includes, upon detecting a condition in the medical infusion pump, displaying a color on a monitor associated to the medical infusion pump.

WO 2008/019015 A2

INTERFACE FOR MEDICAL INFUSION PUMP

This application is being filed on 31 July 2007, as a PCT International Patent application in the name of Smiths Medical MD, Inc., a U.S. national corporation, applicant for the designation of all countries except the US, and
5 William James Evans, a citizen of Great Britain, and Michael L. Blomquist and Rhall E. Pope, citizens of the U.S., applicants for the designation of the US only, and claims priority to U.S. Provisional Patent Application Serial No. 60/835,927, filed August 3, 2006 and U.S. Utility Patent Application Serial No. 11/702,922, filed February 5, 2007.

10

Background

Patients at hospitals and other care centers regularly require controlled drug intake as a part of the patient's prescribed therapy. One form of controlled drug intake is accomplished by infusing fluidic drugs with a medical infusion pump.

15

Medical infusion pumps, in general, provide regulated drug delivery to a patient. These pumps are used to deliver a selected drug or other therapeutic agent to a patient at a predetermined rate that is programmed into the pump. However, programming and managing such pumps can be difficult and cumbersome. Programming typically includes preloading a pump program into a
20 pump and then entering pump parameters or data into the pump through a keypad that is directly in the pump. Each time the pump is programmed, the data must be reentered by hand.

20

Managing the status and locations of pumps also can be difficult. A single pump can be programmed for delivering different fluids in different
25 therapies and in different locations within a hospital. Similarly, the status of a pump and alarms can be difficult to monitor because the pumps are often in locations other than where the caregiver is located and have small displays on which information can be difficult to see.

25

Summary

30

According to a first aspect, an apparatus for indicating a current condition of a medical infusion pump is disclosed. The apparatus includes a memory configured to store two or more color codes, each color code corresponding

to a condition of the medical infusion pump and a color. The apparatus further includes a monitor. The apparatus also includes a programmable circuit in electrical communication with the memory and the monitor, the programmable circuit programmed to generate a graphical user interface displaying a color upon detection
5 of a condition.

According to a second aspect, a method of indicating a current condition of a medical infusion pump is disclosed. The method includes, upon detecting a condition in the medical infusion pump, displaying a color on a monitor associated to the medical infusion pump.

10

Brief Description of the Drawings

FIG. 1 illustrates a pump-computer communication system according to a possible embodiment of the present disclosure;

FIG. 2 illustrates an infusion pump network according to a possible embodiment of the present disclosure;

15

FIG. 3 illustrates the architecture of a computing system that can be used to implement aspects of the present disclosure;

FIG. 4 illustrates the architecture of a pump that can be used to implement aspects of the present disclosure;

20

FIG. 5 is an exemplary infusion pump network according to a possible embodiment of the present disclosure;

FIG. 6A is an exemplary data structure for a pump protocol library according to a possible embodiment of the present disclosure;

FIG. 6B is an exemplary data structure for a pump protocol library according to a possible embodiment of the present disclosure;

25

FIG. 6C is an exemplary data structure for pump protocols according to a possible embodiment of the present disclosure;

FIG. 7 is an exemplary architecture of administrative software for setting global pump protocols according to a possible embodiment of the present disclosure;

30

FIG. 8 is one example of a computer user interface library import screen in accordance with the present disclosure;

FIG. 9 is one example of a computer user interface for administrative software in accordance with the present disclosure;

FIG. 10 is one example of a computer user interface location tab in accordance with the present disclosure;

5 FIG. 11 is one example of a computer user interface therapy tab in accordance with the present disclosure;

FIG. 12 is one example of a computer user interface protocol tab in accordance with the present disclosure;

10 FIG. 13 is one example of a computer user interface drug bar code display screen in accordance with the present disclosure;

FIG. 14 is one example of a computer user interface prescription order form display screen in accordance with the present disclosure;

FIG. 15 is one example of a computer user interface therapy selection screen in accordance with the present disclosure;

15 FIG. 16 is one example of a computer user interface qualifier selection screen in accordance with the present disclosure;

FIG. 17 is one example of a computer user interface drug selection screen in accordance with the present disclosure;

20 FIG. 18 is one example of a computer user interface drug delivery tab in accordance with the present disclosure;

FIG. 19 is one example of a computer user interface weight-based drug delivery tab in accordance with the present disclosure;

FIG. 20 is one example of a computer user interface secondary drug delivery tab in accordance with the present disclosure;

25 FIG. 21 is one example of a computer user interface alarm tab in accordance with the present disclosure;

FIG. 22 is one example of a computer user interface security tab in accordance with the present disclosure;

30 FIG. 23 is one example of a computer user interface appearance tab in accordance with the present disclosure;

FIG. 24 is one example of a computer user interface report tab in accordance with the present disclosure;

FIG. 25 is one example of a computer user interface library export screen in accordance with the present disclosure;

FIG. 26 is a flow diagram of methods and systems for custom programming of a medical infusion pump according to a possible embodiment of the present disclosure;

FIG. 27 is one example of a computer user interface library import screen in accordance with the present disclosure;

FIG. 28 is a flow diagram of methods and systems for editing and loading a protocol for a medical infusion pump according to a possible embodiment of the present disclosure;

FIG. 29 is one example of a computer user interface protocol selection screen in accordance with the present disclosure;

FIG. 30 is one example of a computer user interface therapy selection screen in accordance with the present disclosure;

FIG. 31 is one example of a computer user interface qualifier selection screen in accordance with the present disclosure;

FIG. 32 is one example of a computer user interface drug selection screen in accordance with the present disclosure;

FIG. 33 is a flow diagram of methods and systems for custom programming of a medical infusion pump according to a possible embodiment of the present disclosure;

FIG. 34 is one example of a computer user interface drug delivery customization screen in accordance with the present disclosure;

FIG. 35 is one example of a computer user interface drug delivery customization screen in accordance with the present disclosure;

FIG. 36 is one example of a computer user interface medical infusion pump programming screen in accordance with the present disclosure;

FIG. 37 is a flow diagram of methods and systems for displaying medical infusion pump customizations according to a possible embodiment of the present disclosure;

FIG. 38 is one example of a computer user interface pump comparison screen in accordance with the present disclosure;

FIG. 39 is a schematic front view of a medical infusion pump displaying a change bar according to a possible embodiment of the present disclosure;

FIG. 40 is one example of a computer user interface report generation screen in accordance with the present disclosure; and

FIG. 41 is one example of a computer user interface report screen in accordance with the present disclosure.

Detailed Description

Various embodiments will be described in detail with reference to the drawings, wherein like reference numerals represent like parts and assemblies throughout the several views. Reference to various embodiments does not limit the scope of the claims attached hereto. Additionally, any examples set forth in this specification are not intended to be limiting and merely set forth some of the many possible embodiments for the appended claims.

The following discussion is intended to provide a brief, general description of a suitable computing environment in which the invention may be implemented. Although not required, the invention will be described in the general context of computer-executable instructions being executed by a computer, for example, a hand held computer, a personal computing system, or a medical infusion pump. The structure, creation, and use of a message store hierarchical folder structure are described after the discussion of an exemplary operating environment.

Additionally, the logical operations of the various embodiments of the invention described herein are implemented as: (1) a sequence of computer implemented operations running on a computing system; and/or (2) interconnected machine modules within the computing system. Modules represent functions executed by program code such as commonly available programming languages or as the code found in a dynamic-link library (DLL). The implementation used is a matter of choice dependent on the performance requirements of the pump and the computing systems with which it interfaces. Accordingly, the logical operations making up the embodiments of the invention described herein can be referred to alternatively as operations, modules, and the like.

FIG. 1 illustrates an exemplary embodiment of an infusion pump network 100 having a medical infusion pump 102, a computing system 104, and a communications link 106. The medical infusion pump 102 is configured to deliver therapeutic fluids, such as drugs, saline, or nutrition to a patient. Examples, of
5 medical infusion pumps 102 include ambulatory pumps, stationary pumps, and pole mounted pumps.

The computing system 104 is configured to execute computer-readable instructions, such as computer software. The computing system 104 can be located in a variety of locations such as the point of care (POC) where a patient is
10 being treated, in a healthcare facility at a location remote from the POC, or even at an off-site location remote from the healthcare facility itself. In further embodiments, the medical infusion pump 102 acts as the computing system 104.

In the exemplary embodiment, the computing system 104 is programmed to generate and store pump protocols for execution in the context of a
15 pump application program. Each pump protocol includes a series of pump parameters. Pump parameters refer to settings that define an operational aspect of a medical infusion pump. The pump parameters dictate the control of the pump.

Pump protocols are collections of these pump parameters defining the variable operational characteristics of a medical infusion pump during application of
20 a specific therapy, qualifier, and drug. The pump protocol includes a listing of operational parameters to be included in the pump, and correlates to an index for referring to a specific protocol containing a specific set of pump parameters. The index can be associated with a therapy, qualifier, and drug, and is either contained within the protocol or associated with a specific protocol. The pump protocol
25 includes patient specific pump parameters and non-patient specific pump parameters. Patient specific pump parameters refer to those parameters which are set on a patient-by-patient basis, and for example include the basal delivery rate or bolus amount. Non-patient specific pump parameters refer to those parameters which are set for the pump to perform specific tasks, and do not account for the
30 specific patient to which they are applied. These parameters are generally related to the pump, the infusion pump network, or the medical care to be provided by the pump and/or pump network. Non-patient specific pump parameters can include, for example, a range of permissible values for basal delivery, a range of values and

patterns for basal delivery, a range of permissible values for boluses, a range of values and patterns for extended boluses, a starting value within a particular range of values, alarm values, protocols for data communication, and various flag settings.

5 A pump application program is a program having instructions (e.g., executable code, rules, and/or data) that control operation of the pump for a specific therapy or type of delivery (e.g., continuous delivery, intermittent delivery, pain control, chemotherapy, total parenteral nutrition, etc.). For example, a pump application program might contain instructions that define operation of a pump to accomplish various of the pump parameters. Pump application programs include,
10 for example, pump protocols including both patient specific and non-patient specific pump parameters, and instructions for allocating memory, user interfaces, or algorithms for monitoring various sensors and driving a motor for the pump mechanism.

The communications link 106 connects the pump 102 and computing
15 system 104. In various embodiments, the communications link 106 can include serial or parallel connections, wired or wireless connections, and a direct or networked connection to a computer. Additionally, the pump 102 and the computing system 104 can communicate using any protocol appropriate for data communication. Examples of network connections to a computer include Intranet,
20 Internet, and LAN (e.g., Ethernet). Examples of wired connections to a computer include USB, RS-232, Firewire, and power-line modem connection. Examples of wireless connections include bluetooth, 802.11a/b/g, infrared (IR), and radio frequency (RF).

FIG. 2 illustrates an exemplary embodiment of an infusion pump
25 network 200 having a server 206 networked with a plurality of computing systems 104₁-104_n. The network 200 can be any wired or wireless network that enables data communication between the server, computing systems, and medical infusion pumps. Examples of networks include the Internet, Intranets, and LANs. Each computing system 104 can communicate with a medical infusion pump 102₁-102_n
30 through a communication link 106.

In the exemplary embodiment, the individual computing systems 104_n execute software for generating and managing pump application programs and sets of pump operating parameters. The pump application programs and sets of

pump operating parameters are stored on the server 206 so they can be accessed by other individual computing systems 104_n. The individual computing systems 104_n are also programmed to retrieve previously created pump application programs and sets of pump operating parameters that are stored on the server 206 for viewing,
5 editing, and downloading to medical infusion pumps 102_n.

In alternative embodiments, the medical infusion pumps 102_n can directly access the server to retrieve pump application programs and sets of pump operating parameters. For example, the medical infusion pumps 102_n can be loaded with client software such as a web browser and communicate directly with the
10 network 200, either through a wired or wireless connection as described herein.

In other alternative embodiments, one or more of the computing systems is not configured to communicate directly with a medical infusion pump 102_n, but rather provides administrative access to the server 206 for generating, viewing, and editing pump application programs and sets of pump operating
15 parameters. Additionally, servers, workstations, and other computing systems unaffiliated with the medical infusion pumps 102_n can be included in the network 200.

In yet other alternative embodiments, the software is executed in the server 206. For example, the server functions as an application service provider that
20 communicates user interface and other data entries in mark-up language such as HTML or some other language or protocol that allows a user to execute software from a remote location. In these embodiments, the server 206 can function as an application service provider in which the server provides access to the software for generating and storing pump application programs and pump protocols that a user
25 can create and download to a medical infusion pump. For example, the server 206 could be located at a pump manufacture, pharmaceutical manufacture, pharmacist, or some other third party separate from the user. The server 206 in such an embodiment can be accessed either from an individual computing system 104 or by a medical infusion pump 102 that has networking capabilities and client software.

30 Example embodiments of a server 206 and a medical infusion pump 102 having a web browser are disclosed in United States patent application serial no. 11/066,425, which was filed on February 22, 2005 and is entitled Server for Medical Device, the entire disclosure of which is hereby incorporated by reference.

FIG. 3 illustrates an exemplary architecture that can be used to implement aspects of the present disclosure, including the computing systems 104 and the server 206. The computing system architecture includes a general purpose computing device in the form of a computing system 300. The computing system
5 300 can be used, for example, as the computing system or server of FIG. 2, and can execute program modules included in the administrative software or user software disclosed below.

The computing system 300 including at least one processing system 302. A variety of processing units are available from a variety of manufacturers, for
10 example, Intel or Advanced Micro Devices. The computing system 300 also includes a system memory 304, and a system bus 306 that couples various system components including the system memory 304 to the processing unit 302. The system bus 306 may be any of a number of types of bus structures including a memory bus, or memory controller; a peripheral bus; and a local bus using any of a
15 variety of bus architectures.

The system memory 304 can include read only memory (ROM) 308 and random access memory (RAM) 310. A basic input/output system 312 (BIOS), containing the basic routines that help transfer information between elements within the computing system 300, such as during start up, is typically stored in the ROM
20 308.

The computing system 300 can also include a secondary storage device 313, such as a hard disk drive, for reading from and writing to a hard disk (not shown), and/or a compact flash card 314.

The hard disk drive 313 and compact flash card 314 are connected to
25 the system bus 306 by a hard disk drive interface 320 and a compact flash card interface 322, respectively. The drives and cards and their associated computer readable media provide nonvolatile storage of computer readable instructions, data structures, program modules and other data for the computing system 300.

Although the exemplary environment described herein employs a
30 hard disk drive 313 and a compact flash card 314, other types of computer-readable media, capable of storing data, can be used in the exemplary system. Examples of these other types of computer-readable mediums include magnetic cassettes, flash

memory cards, digital video disks, Bernoulli cartridges, CD ROMS, DVD ROMS, random access memories (RAMs), or read only memories (ROMs).

5 A number of program modules may be stored on the hard disk 313, compact flash card 314, ROM 308, or RAM 310, including an operating system 326, one or more application programs 328, other program modules 330, and program data 332. A user may enter commands and information into the computing system 300 through an input device 334. Examples of input devices might include a keyboard, mouse, microphone, joystick, game pad, satellite dish, scanner, digital camera, touch screen, and a telephone. These and other input devices are often
10 connected to the processing unit 302 through an interface 340 that is coupled to the system bus 306. These input devices also might be connected by any number of interfaces, such as a parallel port, serial port, game port, or a universal serial bus (USB). Wireless communication between input devices and interfaces 340 is possible as well, and can include infrared, bluetooth, 802.11 a/b/g, cellular, or other
15 radio frequency communication systems. A display device 342, such as a monitor or touch screen LCD panel, is also connected to the system bus 306 via an interface, such as a video adapter 344. The display device 342 might be internal or external. In addition to the display device 342, computing systems, in general, typically include other peripheral devices (not shown), such as speakers, printers, and palm
20 devices.

When used in a LAN networking environment, the computing system 300 is connected to the local network through a network interface or adapter 352. When used in a WAN networking environment, such as the Internet, the computing system 300 typically includes a modem 354 or other communications type, such as a
25 direct connection, for establishing communications over the wide area network. The modem 354, which can be internal or external, is connected to the system bus 306 via the interface 340. In a networked environment, program modules depicted relative to the computing system 300, or portions thereof, may be stored in a remote memory storage device. It will be appreciated that the network connections shown
30 are exemplary and other methods of establishing a communications link between the computing systems may be used.

The computing system 300 might also include a recorder 360 connected to the memory 304. The recorder 360 includes a microphone for

receiving sound input and is in communication with the memory 304 for buffering and storing the sound input. The recorder 360 also can include a record button 361 for activating the microphone and communicating the sound input to the memory 304.

5 A computing device, such as computing system 300, typically includes at least some form of computer-readable media. Computer readable media can be any available media that can be accessed by the computing system 300. By way of example, and not limitation, computer-readable media might comprise computer storage media and communication media.

10 Computer storage media includes volatile and nonvolatile, removable and non-removable media implemented in any method or technology for storage of information such as computer readable instructions, data structures, program modules or other data. Computer storage media includes, but is not limited to, RAM, ROM, EEPROM, flash memory or other memory technology, CD-ROM,
15 digital versatile disks (DVD) or other optical storage, magnetic cassettes, magnetic tape, magnetic disk storage or other magnetic storage devices, or any other medium that can be used to store the desired information and that can be accessed by the computing system 300.

 Communication media typically embodies computer-readable
20 instructions, data structures, program modules or other data in a modulated data signal such as a carrier wave or other transport mechanism and includes any information delivery media. The term “modulated data signal” refers to a signal that has one or more of its characteristics set or changed in such a manner as to encode information in the signal. By way of example, and not limitation, communication
25 media includes wired media such as a wired network or direct-wired connection, and wireless media such as acoustic, RF, infrared, and other wireless media. Combinations of any of the above should also be included within the scope of computer-readable media. Computer-readable media may also be referred to as computer program product.

30 FIG. 4 illustrates the architecture of a medical infusion pump 400 that can be used to implement aspects of the present disclosure. A microprocessor 402 is in electrical communication with and controls a pump motor 404, a screen 406, an audible alarm 408, and a vibratory alarm 410. Other embodiments can use a

microcomputer, or any other type of programmable circuit, in place of the microprocessor.

5 The pump motor 404 drives a drive mechanism 412. The drive mechanism 412 delivers the therapeutic fluid to a patient. The drive mechanism can be connected to a plunger system, a peristaltic drive mechanism, or another type of fluid delivery system.

10 The screen 406 can have many different configurations such as an LCD screen. The screen 406 displays a user interface that presents various items of information useful to a patient or caregiver. The audible alarm 408 is a beeper, and an alarm provides actual alarms, warnings, and reminders. Similar to other portable electronic devices such as a cellular telephone, the vibratory alarm 410 provides an alarm to either supplement the audio alarms or replace the audio alarm when an audible beep would be disruptive or not heard. A user can selectively enable or disable the audible 408 and vibratory 410 alarms. In one possible embodiment, however, both the audible 408 and vibratory 410 alarms cannot be disabled at the same time.

20 The microprocessor 402 is in electrical communication with both a random access memory (RAM) 416 and a read only memory (ROM) 418, which are onboard the pump 400 but external to the microprocessor 402 itself. In one possible embodiment, the microprocessor 402 includes internal memory as well. The RAM 416 is a static RAM stores that data that can change over time such as pump settings and a historical log of events experienced by the medical infusion pump 400. The ROM 418 stores code for the operating system and the application programs. The ROM 418 can be any type of programmable ROM such as an EPROM. In one possible embodiment, the RAM 416 has 500 kilobytes of memory capacity and the ROM 418 has 2 megabytes of memory capacity.

30 An infrared (IR) port 420 is in electrical communication with the microprocessor. As explained in more detail below, the IR port 420 provides data communication with an external device such as a computer for programming an application program, programming pump settings, and downloading historical data logs. The medical infusion pump 400 can include other types of communication ports in place of or in addition to the IR port 420. Examples of other possible

communication ports include a radio frequency (RF) port or a port that provides a hard-wired data communication link such as an RS-232 port, a USB port, or the like.

A real-time clock 422 provides a clock signal to the microprocessor 402. An advantage of having a real-time clock 422 is that it provides the program
5 with the actual time in real-time so that the programs executed by the medical infusion pump can track and control the actual time of day that drug delivery and other events occur. Various durations described here are used for alerts, alarms, reminders, and other functions. In one possible embodiment, the timers are formed by the real-time clock 422 and software executed by the microprocessor 402.

10 A keypad 424 also provides input to the microprocessor 402. Although other possible types of keypads are possible, one type of keypad has four buttons and is a membrane-type of keypad, which provides resistance to water and other environmental conditions. The keypad 424 contains soft keys for which the function of the keys can change as a user executes different menu selections and
15 commands.

An audio bolus button 425 optionally provides input to the microprocessor 402. The audio bolus button 425 can program the pump 400 to audibly administer a bolus of drugs or other therapeutic fluids without requiring visual confirmation using the pump. In an example embodiment, the audio bolus
20 button 425 can be pressed a series of times to trigger bolus delivery of a selected volume, based on a preprogrammed trigger granularity. A single button press can represent a bolus of 5 grams, as selected by a user, and subsequent presses of the audio bolus button can represent multiples thereof.

Other inputs into the microprocessor 402 can include an occlusion
25 sensor 426, which is sensitive to occlusions in the therapeutic fluid delivery line; a cartridge sensor 428, which is sensitive to the presence of a therapeutic fluid cartridge; and a motion detector 430, which detects motion of a gear (not shown) in the drive mechanism 412. In an exemplary embodiment, the cartridge sensor 428 includes one or more sensors configured to detect insertion of a therapeutic fluid
30 cartridge. The pump 400 can detect the type of cartridge present via a mechanical interface, and can include in the pump software instructions regarding operation in conjunction with the cartridge. Examples of cassette sensing features are described, for example, in United States Patent No. 5,531,697, filed on April 15, 1994, issued

on July 2, 1996, and entitled Systems and Methods for Cassette Identification for Drug Pumps.

FIG. 5 illustrates a schematic architecture of a medical infusion pump network 500 according to an exemplary embodiment. The medical infusion pump network 500 includes an administrator computer 502 communicatively connected to the server 206 of FIG. 2, which includes a database 504. The medical infusion pump network 500 also includes one or more medical infusion pumps 102 and computing systems 104.

The administrator computer 502 and computing systems 104 are systems such as those described above in conjunction with FIG. 3. The administrator computer 502 includes administrative software installed on or accessible to the computer for generating one or more libraries 508 of pump protocols 510. An exemplary embodiment of the administrative software is described below in FIGS. 7-25.

In the present disclosure, libraries refer to collections of pump protocols generated using the administrative software described herein. Libraries can be stored in files, databases, or other data structures. Libraries contain pump protocols as well as indices pointing to the protocols, and are loaded in user software to select a specific pump protocol for operation of a medical infusion pump.

The computing systems 104 include user software for accessing one or more libraries 508 of protocols 510 and programming a medical infusion pump 102 with a protocol 510 or a library 508. In one possible embodiment, the computing systems 104 are optional in that the user software resides directly on the medical infusion pumps 102. An exemplary embodiment of the user software is described below in FIGS. 26-41.

The medical infusion pumps 102 connect either to a computing system 104 or directly to the server 206, and are described above in conjunction with FIG. 4. In a first embodiment, the medical infusion pumps 102 are configured to accept a pump protocol from the server 206 or the computing system 104. In a second embodiment, the medical infusion pumps 102 are configured to accept a library 508 of pump protocols 510 directly from the server 206 or from the computing system 104.

The database 504 contains pump protocol data 506 and log files 516. The pump protocol data 506 forms a plurality of libraries 508 which in turn each include a number of protocols 510. Each protocol 510 is stored as a data record, and includes a set of parameters, including patient specific pump parameters 512a and
5 non-patient specific pump parameters 512b, as described above. Each library 508 can contain one or more pump protocols 510.

The log files 516 include log data regarding access and usage of the libraries 508, and can include information related to the administrator computer 502, the medical infusion pumps 102, or the computing systems 104 authorized to
10 connect to the server 206. In one possible embodiment, the log files include access records, which record instances in which medical infusion pumps access a library 508 on the server 206.

FIGS. 6A-6C illustrate exemplary data sets accepted by medical infusion pumps 102. The specific data set accepted by a medical infusion pump 102
15 is dependent upon that pump, but in various embodiments, the pumps accept data sets representing pump libraries, pump protocols, or pump programs. FIG. 6A shows a library 508 in greater detail. The library 508 can be loaded into the memory of a medical infusion pump 102, allowing a user of the pump to select a protocol for operation. The library 508 includes a number of protocols 510, which include
20 patient specific pump parameters 512a, non-patient specific pump parameters 512b, and an index 514. The number of protocols 510 within a given library can vary, and depends upon the number defined using the administrative software.

The total number of pump parameters 512a-512b remains constant for each particular model of pump, but can vary between types of pumps.
25 Additionally, the pump parameters 512a-512b can be configured in a number of formats within each protocol 510 for the same type of pump. For example, the number of patient specific pump parameters 512a can vary between protocols due to the specific type of drug and therapy applied. For example, a protocol defining a continuous drug delivery may only require a single patient specific protocol, namely,
30 the drug delivery rate. In another example, a protocol defining intermittent drug delivery may require additional patient specific pump parameters, such as the time between drug delivery phases, a bolus amount, patient bolus amounts, and other parameters. The number of non-patient specific pump parameters 512b represents

the difference between the total number of parameters programmable into a pump and the number of patient specific pump parameters 512a as dictated by the therapy and drug applied.

The index 514 can be any generic index referencing a specific location within the library. Each index is unique within the library, although another library may contain the same index and relate that index to a different set of pump parameters contained within that library. In the exemplary embodiment shown, the index 514 includes therapy, qualifier, and drug regions. By selecting a combination of a therapy, a qualifier, and a drug, a user of the system can select one of the protocols 510 from the library 508. Therapies, as referred to herein, are the methods of patient treatment for diseases or generalized rehabilitation. For example, a therapy can be an epidural treatment or patient-controlled analgesia. Qualifiers include factors affecting the administration of a therapy, such as weight, age, or sensitivity of a patient to a specific therapy. Drugs refer to any therapeutic fluids deliverable by a medical infusion pump.

Each of the protocol entries on the server can be assigned an identification code in order to ensure that the medical infusion pumps access a correct library and/or protocol, and that the protocols on the medical infusion pumps and computing systems associated with the server 206 of FIG. 5 are up to date. These identification codes are generated by the server and stored in conjunction with the protocol in the server 206, as well as in the medical infusion pump 102 and/or its associated computing system 104. The identification codes can be generated using globally unique identifiers (GUID), and are used to track the specific protocol and/or library accessed by each pump 102 or computing system 104 in the database 504. A GUID is a 128-bit pseudo random number used to provide a statistically unique identifier for corresponding the protocol on the pump 102 to the protocol as stored in the server 206. The GUID can be generated by the server 206 and transmitted alongside the protocol and/or library when transmitted to the computing system 104 or infusion pump 102. The server copy of the GUID can be stored alongside the protocol in the database 504, and the pump copy can be stored in RAM in the medical infusion pump 102 or computing system 104. Each time the medical infusion pump 102 or computing system 104 accesses the protocol on the server 206, the GUID assigned to that protocol as stored in the pump 102 or computing

system 104 can be matched to the protocol as stored in the database 504 to ensure that the correct protocol is accessed. Different protocols in different libraries can use the same index of therapy, qualifier, and drug, and can be stored in the same database, but have a unique GUID and are therefore uniquely identifiable. In a possible embodiment, the GUID system can be used in conjunction with user access control systems, such as are disclosed in conjunction with both the user software and administrative software, below. In a further embodiment, a new GUID can be generated and associated with each protocol when first created, or optionally each time the protocol is edited using administrative software, such as is described below.

FIG. 6B illustrates a second possible embodiment of a library 508. The library 508 contains a number of protocols 510, which, as in FIG. 6A, include patient specific pump parameters 512a, non-patient specific pump parameters 512b, and an index 514. The number of protocols 510 included in each library 508 may vary, so the number of protocols 510 in the library 508 of FIG. 6B may be different from the number of protocols in the library shown in 6A. In the embodiment shown, the library 508 contains one protocol 510. Additionally, one or more of the therapy, qualifier, and drug regions in the index 514 can be replaced by other index criteria, such as locations, pump programs, doctor identification, or other indexable criteria capable of referring to a unique pump protocol within the library. In the example shown, the "therapy" region is used to select a pump program, such as a continuous delivery program, an intermittent delivery program, or a specific type of program such as a pain management program. The "qualifier" region is used to select a doctor, and may be the name of a doctor using the infusion pump network of FIG. 5.

FIG. 6C illustrates a series of possible pump protocols 510a-510c. The pump protocols 510a-510c incorporate patient specific pump parameters 512a and non-patient specific pump parameters 512b. Protocols 510 are specific to a pump 102, in that the pump has a specific number and type of parameters that are programmable. Therefore, the total number of pump parameters remains constant. However, the number of patient specific pump parameters 512a can vary depending upon the protocol selected for programming into the pump, which in turn dictates that the number of non-patient specific pump parameters 512b varies as well. In a possible embodiment, one of the protocols 510 is selected using a computing system

104 or infusion pump 102. If selected using the computing system, the protocol is then programmed into the medical infusion pump 102.

In another possible embodiment, the pump protocol 510 is selected using the infusion pump 102 or the computing system 104. The pump protocol 510 is then incorporated into a pump program to provide a set of instructions dictating the operation of a medical infusion pump 102 according to the protocol 510. The complete pump program is then downloaded into the pump 102. In yet another possible embodiment, the pump program is downloaded to the pump 102 at a different time from the pump protocol 510. In still a further embodiment, multiple pump programs reside within the pump 102, and the pump protocol 510 contains a parameter which dictates which pump program is to be used. In a further embodiment, the pump program within the pump is altered based on one or more of the pump parameters included in the pump protocol 510.

FIG. 7 illustrates exemplary architecture of administrative software 700 for generating one or more libraries of pump protocols. The software 700 can operate within the server 206, pump 102, computing system 104, or a combination thereof.

The administrative software 700 allows a user, for example a doctor, nurse, pharmacist, or other caregiver, to create, define, and edit pump application programs and protocols for execution in and control of medical infusion pumps 102. For example, the administrative software 700 can generate protocols and programs that can be loaded using the user software described in FIGS. 26-41, below.

The administrative software 700 provides protocol-based programming of medical infusion pumps in which the user creates a pump application program by designating a particular therapy and other criteria such as a location and qualifiers (e.g., patient age, weight, skin surface area). Once criteria are selected, the administrative software 700 applies rules and other logic that assembles sets of pump parameters into a pump protocol. For example, the administrative software 700 might be used to select one delivery pattern and enable bolus delivery if the selected therapy is for delivering pain medication and another delivery pattern and not enable bolus delivery if the selected therapy is for parenteral nutrition. In another example, the administrative software 700 might be used to select one range of permissible delivery rates if one of the criteria indicates the patient is an

adolescent and different range of permissible delivery rates if the patient is an adult. Other embodiments permit programming a medical infusion pump 102 without using therapy-based programming. Additional embodiments of protocol- or therapy-based programming is discussed in more detail in United States patent
5 application serial no. 11/003,147, filed on December 3, 2004 and entitled Programming Medical Pumps with Electronic Standing Order Template, the entire disclosure of which is hereby incorporated by reference.

Operation of the software 700 begins at a start module 702. The start module 702 corresponds to initial execution of the administrative software by
10 clicking on an icon on the computer or by some other mechanism for executing software. Upon startup, the software 700 connects to a library loaded in the database 504 of FIG. 5.

Following the start module, operational flow optionally proceeds to a load library module 704, which allows a user to access a listing of library files
15 available to the administrative software 700. The library files contain one or more libraries, which in turn contain a collection of pump protocols as described above. The collection of library files can be stored in the server 206 or in one or more individual computing systems 102. The load library module 704 allows the user to select a library file containing one or more libraries for viewing, editing, and
20 downloading to a medical infusion pump 102. If a user does not want to download or otherwise access an existing library, it can selectively bypass the load library module 704. An example of when a user bypasses the load library module 704 is if the user plans to only create a new library or edit one or more protocols within the currently loaded library. In an alternative embodiment, the software always
25 executes the load library module 704 and the user then selectively chooses whether to load any previously created libraries via a stored library file.

Following the start module 702 and optional load library module 704, operational flow proceeds to a login module 706. The login module 706 regulates user rights in the software 700. User rights define access levels to the currently
30 connected library in user software, and are configurable for users such as doctors, nurses, or other caregivers. Based on the user rights assigned to a caregiver, that user will have a set access level allowing the user to view, add, or edit pump libraries within the user software, described in detail below. Access levels can be

set according to a variety of criteria. Examples include the type of caregiver (e.g., physician, nurse, pharmacist), location (e.g., hospital, clinic, pharmacy, manufacturer), or a particular department within a location.

In possible embodiments, different access levels also provide
5 different rights with respect to a particular pump protocol or pump operational parameters. For example, one access level might give a user a right to edit, create, and download pump protocols and/or pump application programs. One access level might permit a user the right to edit, create, and download only specified pump
10 parameters, such as the patient specific pump parameters described in conjunction with FIGS. 5-6. One access level might permit a user the right to only edit or download pump parameters. One access level might permit a user the right to only view and download pump parameters. Different embodiments can include the ability to provide an access level for a user any combination of rights to create, edit,
15 parameters. An example of lock levels are disclosed in United States Patent 6,475,180, issued on November 5, 2002 and entitled Drug Pump Systems and Methods, the entire disclosure of which is hereby incorporated by reference.

Once the user is logged in and the library is optionally loaded, the user selectively executes three different modules, a library module 708, a therapy
20 module 710, and a protocol module 712.

The library module 708 assigns a label that identifies an entity and user attributes for the selected entity. Entity attributes can be properties specific to the library, such as a name of a doctor, a name of a healthcare regimen, or a location of the medical infusion pump or pump network, for example the hospital or
25 department at which the pump is located. User attributes define the users allowed to access and modify pump parameters for protocols associated with a particular library by using the medical infusion pump, and can also define users allowed to modify pump protocols using user software, as described below. The library module 708 contains a library definition module 714 and a user rights module 716, which are
30 configured to perform these tasks, respectively.

The therapy module 710 adds and modifies therapies, qualifiers associated with the therapies, and drugs. The therapy module includes a therapy definition module 718, a qualifier definition module 720, and a drug definition

module 722. The therapy definition module 718 controls addition and editing of therapies, which are the methods of patient treatment for diseases or generalized rehabilitation as previously described. The qualifier module 720 defines qualifiers and associates the qualifiers with one or more therapies. The drug definition module
5 722 defines one or more drugs that can be used in the medical infusion pump.

The protocol module 712 adds, edits, and defines protocols by associating therapies, qualifiers, and drugs with pump parameters to form libraries of pump protocols for a medical infusion pump. The protocol module 712 allows a user to select a therapy defined in the therapy definition module 718. The protocol
10 module further allows the user to associate a qualifier defined in the qualifier definition module 720 with the selected therapy. For example, one or both of "adults" and "children 5-10 years" qualifiers can be associated with an epidural therapy. The protocol module 712 also associates one or more drugs with each therapy and qualifier combination, indicating that use of the drug is appropriate for
15 that therapy and qualifier. The protocol module 712 guides the user in defining a protocol by assigning default pump parameters to be associated with the selected therapy, qualifier, and drug.

The protocol module 712 allows a user to associate more than one qualifier to each therapy, and also allows a user to associate more than one drug to
20 each therapy and qualifier combination. For example, a protocol used in an epidural therapy for an adult can include a higher basal delivery rate parameter than a protocol used with a child for the same therapy. Likewise, usage of one drug can require a higher or lower dosage than another drug for the same therapy and qualifier because of concentration, reaction, or other factors.

Operational flow proceeds from the modules 708, 710, 712 to an optional export library module 724. The export library module 724 saves the defined or edited pump application programs and parameters in a file or other data
25 structure that can be loaded by the administrative software 700 at another time or location, or can be loaded by user software such as described below in FIGS. 26-41. If the export library module 724 does not execute, the library remains within the
30 normally connected database 504 of FIG. 5, but is not extracted into a file for portability to a pump not connected to the network or to another medical infusion pump network. One or more libraries can be exported into a single library file.

Operation of the software 700 terminates at an end module 726. The end module 726 corresponds to termination of the administrative software 700 by clicking on a close window button on the computer or by some other mechanism for terminating execution of software.

5 In one embodiment of the administrative software 700, a user of the software 700 defines each protocol included in a library. In defining each protocol, the user assigns the index to the protocol, such as the therapy, qualifier, and drug defined in the modules 708, 710, 712 above. In a second possible embodiment, the administrative software includes a number of default settings or pump parameter
10 modifications used when specific therapies, qualifiers, or drugs are selected. The user selects a therapy, qualifier, and drug to associate with a pump protocol. The administrative software 700 can include instructions dictating that selection of one or more of the therapies, qualifiers, and drugs sets or modifies one or more of the patient specific pump parameters or non-patient specific pump parameters. In one
15 example of this second embodiment, a user setting a drug having a maximum safe consumption rate will trigger the administrative software 700 to preset an acceptable range of programmable delivery rates and a default delivery rate in the protocol, as well as alarms or other non-patient specific pump parameters. In another possible
20 example of this second embodiment, a user setting a qualifier indicating a low age, such as "Children 5-10 years old", will set or adjust the protocol to result in a low delivery rate and demand dose being incorporated into the protocol, and will set one or more parameters related to alarms for use in a medical infusion pump.

Referring now to FIG. 8, a library import screen 800 is shown. The library import screen corresponds to the load library module 704, and is optionally
25 used to load a library file containing one or more libraries of pump programs in the administrative software 700 of FIG. 7, above. The library import screen 800 includes a file selection field 802 and a password field 804. The term field as used herein can include the window, or screen, generally, and can also include menus, selectable lists, or buttons within the window.

30 The file selection field 802 displays one or more library files that are available to be selected. The file selection field 802 allows a user to select one or more of the protocol library files, in conjunction with selection control buttons 806a, 806b. The password field 804 controls access to the selected library file by

requiring a user to input a correct password associated with the selected file, or library within the file.

A location field 808 presents the directory path of the selected library file and a browse button 810 provides browsing capabilities to allow a user to find a library file other than those displayed in the library selection field 802.

In use, the library import screen 800 initially presents a listing of library files in the file selection field 802 available to the administrative software. The user selects one of the displayed library files, or presses the browse button 810 to search for additional library files. The user selects a library file by clicking on the displayed library file, or by other selection method. The directory path of the selected library file appears in the location field 808, and the user then enters a password in the password field 804 corresponding to the selected library file. The user confirms the choice using the selection control buttons 806a, 806b. The administrative software confirms that the user-entered password is correct, and accesses the library file.

Referring now to FIG. 9, a user interface 900 for the administrative software 700 of FIG. 7 is shown. The user interface 900 includes features associated with the modules incorporated in the administrative software 700, and corresponds to the screen displayed following the start operation 702 of FIG. 7. The user interface 900 includes a number of features providing global control and status information to a user.

Global control features included in the user interface 900 relate to library and pump access settings. A library field 902 provides a listing of libraries currently loaded by the software. The library field 902 contains the libraries which have been loaded using the load library module 704 of FIG. 7. A login button 904 generates a login screen that checks whether a user has the right to modify pump parameters and protocols.

The user interface 900 can display information related to the status of the network of medical infusion pumps. A database identification field 906 identifies the database 504 currently connected to the administrator computer 502, as shown in FIG. 5. In additional exemplary embodiments, additional information can be displayed, such as the location of the server and/or medical infusion pumps

connected to the infusion pump network or a serial number or other identification number of the pumps or libraries.

The user interface 900 also includes a location tab 920, a therapy tab 940, and a protocol tab 960. Referring back to FIG. 5, the location tab 920
5 corresponds to the location module 708, the therapy tab 940 corresponds to the therapy module 710, and the protocol tab 960 corresponds to the protocol module 712. By clicking on or otherwise selecting a tab, a user transfers operation to the module corresponding to that tab.

Referring now to FIGS. 10-12, the user interface 900 is shown, and
10 location tab 920, therapy tab 940, and protocol tab 960 are described in greater detail. FIG. 10 shows the user interface 900 after selection of the location tab 920, or after initial execution of the login module 706 and optional load library module 704.

A region of the location tab 920 supersedes the therapy tab 940 and
15 protocol tab 960. Within the region exposed by the location tab 920, a library selected in the library field 902 populates a library description field 1002 and a user accounts field 1004. The library description field 1002 describes the library that is currently loaded, and includes attributes of the location in which the library is used, or other information about the library. The location attributes can include the name
20 of the hospital or the department in the hospital associated with the library. The additional information can include information related to the users of the library, the contents of the library, or other information. The library description field 1002 corresponds to the location attributes module 714 of FIG. 7. The user accounts field 1004 displays a listing of user accounts associated with the library. Each of the user
25 accounts generally correspond to a user, but can also correspond to a location of a computing system 104 or medical infusion pump 102. The user accounts define the persons allowed to access and modify pump protocols, and/or parameters associated with the library at the point of care, such as a medical infusion pump 102 or computing system 104 executing user software. The user accounts field 1004
30 enables a user to add, edit, or delete users from the listing using the buttons 1006a-1006c provided. The user accounts field corresponds to the user rights module 716 of FIG. 7.

In a possible embodiment, a log report field 1008 can optionally direct the server 206 to generate a record for various events occurring in the database 504. For example, a log report can be created each time a protocol is sent to a medical infusion pump 102. Alternately, a log report is created each time a library is sent to a medical infusion pump 102. In further embodiments, a log report is created each time a library is edited or accessed. In still further embodiments, globally unique identifiers are entered into the log report related to instances where computing systems 104 and/or infusion pumps 102 access a library in the database 504. In the embodiment shown, the log report field 1008 is a selectable check box, but can be implemented as any other type of selectable field.

A password field 1010 sets a password for the currently selected library file. The password protects access to the library from the perspective of user software residing on either a computing system 104 or an infusion pump 102 such that only users with knowledge of the password associated with the library file can load the library using the load library module 704 of FIG. 7 and user interface 800.

Referring now to FIG. 11, the user interface 900 is shown with the therapy tab 940 selected. A region of the therapy tab 940 supersedes the location tab 920 and the protocol tab 960. The therapy tab 940 corresponds to the therapies module 710, which defines therapies, qualifiers, and drugs used by medical infusion pumps associated with the administrative software 700 of FIG. 7. The therapy tab 940 includes a therapy definition field 1102, a qualifier definition field 1104, and a drug definition field 1106.

The therapy definition field 1102 corresponds to the therapy definition module 718 of FIG. 7, and includes a therapy listing 1108, a therapy notes field 1110, and control buttons 1112a-1112c. The therapy listing 1108 displays the therapies currently defined in the library displayed in the library listing 802. Two example therapies, "Patient Controlled Analgesia" and "Epidural", are shown in the therapy listing 1108, and are associated with a library named "NeoNatal Intensive Care Unit". The therapy notes field 1110 contains administrative user-defined notes related to the therapy selected in the therapy listing 1108. The notes relate to a therapy, and include information related to administration of the therapy, such as messages related to dosage, bolus amount, or administration. The therapy notes field 1110 presents administrative user-created notes to convey therapy-specific

information to caregivers using or programming the pump. Control buttons 1112a-1112c allow a user to add, edit, and delete therapies from the therapy listing 1108 and therapy notes field 1110.

The qualifier definition field 1104 corresponds to the qualifier definition module 720 of FIG. 7, and includes a qualifier listing 1114, a qualifier notes field 1116, and control buttons 1118a-1118c. The qualifier definition field 1104 lists the qualifiers associated with the therapy selected from the list shown in the therapy listing 1108. For example, an administrative user might add a "Children 5-10 yrs" entry to the qualifier listing 1114. Each therapy relates to one or more qualifiers, such as a general qualifier, a weight based qualifier, or an age based qualifier. The qualifier notes field 1116 contains notes describing the difference in application of the therapy based on the qualifier selected. In the case of the "Children 5-10 years" qualifier, the qualifier notes field 1116 can indicate a lower than normal dosage to be administered. Notes associated with the qualifier display in the qualifier notes field 1116 only when the qualifier is selected.

The drug definition field 1106 corresponds to the drug definition module 722 of FIG. 7, and includes a drug listing 1120 and control buttons 1122a-1122d. The drug listing 1120 includes information related to the therapeutic fluid used in the medical infusion pump. This information can include the drug name, identification code, units, concentration, and usage. The control buttons 1122a-1122c allow a user to add, edit, or delete drugs in the drug listing 1120. A bar code control button 1122d generates a bar code screen including information related to a drug.

Referring now to FIG. 12, the user interface 900 is shown with protocol tab 960 selected. The protocol tab 960 supersedes the location tab 920 and the therapy tab 940. The protocol tab 960 corresponds to the protocol module 712 of FIG. 7, and associates therapies, qualifiers, and drugs to allow a user to define protocols by setting pump parameters for inclusion within pump programs. The protocol tab 960 includes a protocol field 1202 and control buttons 1204a-1204e. The protocol field 1202 lists the combinations of therapies, qualifiers, and drugs for which a protocol is defined. The protocol field 1202 also lists a pump type 1203 for which each protocol is defined. By specifying a pump type for which each protocol is defined, it is possible to enable pump-specific protocol programming, while still

using the administrative software 700 and user software, described below, for all pump types configurable using the protocol-based programming scheme described herein.

5 The control buttons 1204a-1204e operate to add, edit, view, and/or delete the protocol for the combinations of therapies, qualifiers, and drugs. The control buttons 1204a-1204d allow a user to set pump parameters so as to define the protocol. Control button 1204e generates a prescription form screen representing the protocol for the selected therapy, qualifier, and drug.

10 FIG. 13 shows a bar code screen 1300 displaying a bar code for a drug defined in the administrative software 700 of FIG. 7. A drug identification code is associated with each drug in the drug listing 1120 of FIG. 11. A drug is selected from among those in the drug listing 1120 of FIG. 11, and the bar code screen 1300 displays upon clicking on the control button 1122d.

15 The bar code screen 1300 includes a print preview field 1302, a printer drop down menu 1304, a label information menu 1306, a copies drop down menu 1308, and control buttons 1310a-1310c. The print preview field 1302 displays a bar code associated with the selected drug. The printer drop down menu 1304 lists available printers configured to print the bar code. The label information drop down menu 1306 defines the label configuration to which the bar code is directed. The
20 label configuration includes, for example, the size and layout of the label paper. The copies drop down menu 1308 dictates the number of copies of the bar code that are printed. Control buttons 1310a-1310c provide printing, cancellation, and help procedures to a user.

25 In use, the bar code corresponds to a drug identification code associated with the drug. In a possible embodiment, a pharmacist can print the barcode associated with the drug using the administrative software 700 via the bar code screen 1300 and affix the printed bar code label to the drug container. The labels indicate the drug and dosage being delivered by the pump. This provides an easily accessible and visually prominent indication of the drug being delivered by a
30 medical infusion pump.

A caregiver connecting a drug to a medical infusion pump 102 will scan the bar code on the drug container, which will correspond to the drug identification code associated with the drug in the administrative software 700. This

ensures that the caregiver affixes the drug to the pump which corresponds to the protocol selected using administrative software 700.

FIG. 14 shows a prescription form screen 1400 displaying prescription information related to a specific therapy, qualifier, and drug. The prescription form screen displays upon selection of control button 1204e on protocol tab 960. The prescription form screen 1400 corresponds to the currently selected protocol, and therefore incorporates information specific to the selected protocol. Prescription information includes, for example, directions for application of the drug according to the defined protocol, and parameters such as drug delivery rates or bolus amounts. The prescription form also includes usage notes describing operation or application of the selected therapy, qualifier, and drug.

A doctor using the prescription form dictates the protocol used in the pump by using a prescription form analogous to the prescription form screen. The prescription form 1400 generated by the administrative software will correspond to the prescription form completed by the doctor. In the embodiment shown, the prescription form screen 1400 that is generated includes information specific to the drug and therapy administered, which may be notes related to administration of the drug and therapy as dictated by the doctor. For example, the prescription form will include drug information, such as the name, type, concentration, and notes regarding the drug, and will also include information related to patient specific pump parameters associated with a selected therapy, qualifier, and drug selected from a library.

Referring now to FIG. 15-17, a protocol definition user interface 1500 defines the relationships between the therapies, qualifiers, and drugs added to the library by the therapy module 710 and using the therapy tab 960. A user defines a protocol by setting an index by sequentially selecting a therapy, a qualifier, and a drug, and then by associating pump parameters with that index. In setting the index using the user interface 1500, a therapy drop down menu 1502 connects to a therapy notes field 1504, a qualifier drop down menu 1506 connects to a qualifier notes field 1508, and a drug drop down menu 1510 connects to a drug notes field 1512.

FIG. 15 shows the initial state of the protocol definition user interface 1500. A list of therapies appears in the therapy drop down menu 1502, while the remaining menus 1506, 1510 and fields 1504, 1508, 1512 remain blank.

FIG. 16 shows the state of the protocol definition user interface 1500 after a therapy is selected in the therapy drop down menu 1502. Notes related to the selected therapy appear in the therapy notes field 1504, and a list of qualifiers associated with the selected therapy appears in the qualifier drop down menu 1506.

- 5 The notes correspond to the therapy notes entered in the therapy notes field 1110. The list of qualifiers corresponds to the qualifiers associated to the therapy by listing in the qualifier listing 1114. For example, two qualifiers, "Adults & Children over 5 yrs" and "Children 3-5 yrs" can be associated with the "Patient Controlled Analgesia" therapy. When that therapy is selected in the therapy drop down menu
- 10 1502, the two associated qualifiers populate the qualifier drop down menu 1506.

- FIG. 17 shows the state of the protocol definition interface 1500 after a qualifier is selected in the qualifier drop down menu 1506. Notes related to the qualifier appear in the qualifier notes field 1508, and a list of drugs available within the library or database appears in the drug drop down menu 1510. The notes
- 15 correspond to changes in the therapy due to the qualifier selected. Continuing the example from FIG. 14, the "Adults & Children over 5 yrs" qualifier is selected. Notes related to customization of the "Patient Controlled Analgesia" therapy based on the selected qualifier are displayed in the qualifier notes field 1508. Furthermore, a list of four drugs/drug concentrations appear in the drug drop down menu 1510,
- 20 which are the drugs defined and available within the library or database.

Analogously, upon selection of a drug from the drug drop down menu 1510, drug notes (not shown) appear in the drug notes field 1512.

- Referring now to FIGS. 18-24, a parameter user interface 1800 provides additional tabbed screens allowing a user to define pump parameters. The
- 25 pump parameters complete the protocol definition associated with the therapy, qualifier, and drug combination selected in FIGS. 15-17. As described in conjunction with FIG. 7, selection of the index defined by a therapy, qualifier, and drug using the protocol definition interface 1500 may trigger the administrative software 700 to set or modify one or more of the pump parameters. Alternately,
- 30 none of the pump parameters are set during the definition and selection of the index formed by the therapy, qualifier, and drug using the protocol definition interface 1500. In such an embodiment, the parameter user interface 1800 is used to define the pump parameters that are programmable into a medical infusion pump.

The parameter user interface 1800 includes a status region 1802, a protocol activation field 1804, and control buttons 1806a-1806c. The status region 1802 displays the therapy, qualifier, and drug associated with the assignable pump parameters. The protocol activation field 1804 publishes the protocol within the library such that the protocol is visible to user software accessing the library when it resides within the database 504 of FIG. 5. The protocol activation field 1804 allows a user of the administrative software 700 to control when protocols become available for use by user software. In some circumstances, it can be advantageous to prevent user software from accessing data, particularly while that data is being edited in the administrative software 700. One example of user software used to access pump protocols is illustrated below in conjunction with FIGS. 26-41. The control buttons 1806a-1806c provide save, cancel, and help options to a user of the administrative software.

The parameter user interface 1800 further includes a number of tabs, including a drug delivery tab 1810, a secondary drug delivery tab 1820, an alarm tab 1830, a security tab 1840, a display/sound tab 1850, and a report tab 1860. Parameters set within each of the tabs are discussed in FIGS. 18-24.

FIG. 18 shows the parameter user interface 1800 with the drug delivery tab 1810 selected. The drug delivery tab 1810 allows a user to set drug delivery rate parameters, and includes a general settings region 1812 and a patient specific parameters region 1814.

The general settings region 1812 includes verification settings 1816 and weight based settings 1818. The verification settings 1816 includes drug verification and caregiver verification settings. Specifically, the verification settings require that a caregiver verifies that the correct drug is provided to the medical infusion pump. The verification settings also require a second caregiver to verify the settings of the medical infusion pump. The weight based settings 1818 set a weight based protocol at a programmable, variable weight limit. By weight based protocol, it is intended that dosage delivery rates, boluses, thresholds, and other delivery parameters change from a "dosage per hour" basis to a "dosage per weight factor" rate, where the weight factor can be on a per unit measure weight basis for the user of the medical infusion pump 102 (i.e. "per kilogram" or other), or based on the user's body surface area, a weight based therapy, or other options.

The patient specific parameters region 1814 includes a continuous rate region 1822, a demand dose region 1824, and a demand dose lockout region 1826. Continuous rate refers to the constant drug delivery rate of the medical infusion pump, also referred to as the basal rate. Demand dose refers to an added drug delivery bolus amount delivered by the pump upon a demand by a patient. Demand dose lockout refers to the time interval after a demand dose is delivered, during which another demand dose will not be delivered by the pump.

The continuous rate region 1822 includes a meter, shown as a slider bar 1828 and an indicator 1829. The meter generally has two or more locations, each corresponding to a parameter value that can be programmed in the medical infusion pump. Generally, the positional relationship of the meter indicates the setting of the meter. In a possible embodiment of the slider bar 1828 shown, the indicator 1829 is movable relative to the slider bar 1828 to set a default value, or "initial value" continuous drug delivery rate parameter. In a second possible embodiment, the default value is set using an initial value gauge 1832.

The continuous rate region also includes hard limit gauges 1834, soft limit gauges 1836, and user interface options 1838a-1838c. The initial value gauge 1832, hard limit gauges 1834, and soft limit gauges 1836 include values, which may include numerical ranges. The hard limit gauges 1834 set a hard maximum and hard minimum which form an acceptable pump programming range. The range of acceptable pump activity represents the absolute maximum and minimum values programmable into the pump by user software, as described below. This configuration allows for control of the range of values visible to a user of the medical infusion pump or associated computing system.

The limits set by the soft limit gauges 1836 represent a manually exceedable threshold value. The soft limit can be overridden by a user of a medical infusion pump on a pump-by-pump basis. Pump activity outside the range defined by soft limits can trigger an alarm or otherwise alert a caregiver that a pump is functioning outside of the usual operational range of the pump. A variety of alarm levels or alerts can be set by the soft limit gauges 1836. For example, the alert can be a flag set in the software. The alert could additionally be an audible alarm, or a visual indicator displayed on at least a portion of the medical infusion pump. The

visual indicator could be a flashing indicator or changed/changing color on the display of the medical infusion pump.

In a second possible embodiment, the hard limit gauges set a non-limiting range, and the user software described below can be programmed within its full operational range. In such an embodiment, pump activity outside the range set by the hard limit gauges 1834 can trigger an alarm or otherwise alert a caregiver that a pump is functioning outside of the usual operational range of the pump. In this embodiment, the soft limit gauges 1836 set a narrower range, operation outside of which can trigger a warning or second alarm indicating pump activity outside of an expected range of pump operation. This warning or second alarm indicates a pump condition less serious than the alarm triggered by the hard limits.

The user interface options 1838a-1838c enable the user software to display and edit the delivery rate, editing of the delivery rate, and requiring comments by users of a pump who wish to exceed the soft limits when setting the drug delivery rate. Selection of the display option 1838a publishes the pump parameter so that the value is visible to a user of the pump or computing system.

The demand dose region 1824 includes a slider bar 1842 and an indicator 1843. The slider bar and indicator operate in a similar manner to the slider bar 1828 and indicator 1829 in the continuous rate region 1822, but control demand dose settings. Likewise, the demand dose region 1824 includes an initial value gauge 1844, as well as hard limit gauges 1846 and soft limit gauges 1848 setting visible thresholds and triggering alarms as in the continuous rate region 1822. Demand dose options 1852a-1852c provide analogous display, editing, and comment options to the user interface options 1838a-1838c.

The demand dose lockout region 1826 includes a slider bar 1854 and an indicator 1855, and also includes an initial value gauge 1856, hard limit gauges 1858, and soft limit gauges 1862. Each of these features functions analogously to those discussed above in conjunction with the continuous rate region 1822. The demand dose lockout region also includes lockout options 1864a-1864c analogous to the user interface options 1838a-1838c.

FIG. 19 shows the parameter user interface 1800 with the drug delivery tab 1810 modified to provide a weight based drug delivery protocol. The modification of the user interface 1800 occurs in the drug delivery tab 1810 upon

user selection of the weight based settings 1818 discussed in FIG. 18. The continuous rate region 1822 and demand dose region 1824 are modified to reflect dosage rates on a "per kilogram" or other weight measure basis.

FIG. 20 shows the parameter user interface 1800 with the secondary drug delivery tab 1820 selected. The secondary drug delivery tab 1820 provides additional medical infusion pump programming options for assigning pump parameters in a specific protocol. The secondary drug delivery tab 1820 includes a dosing limit region 2002, a patient specific parameter region 2004, a titration region 2006, a maximum delivery rate gauge 2008, and a maximum clinician bolus gauge 2010.

The dosing limit region 2002 displays limits for total drug delivery within a specified amount of time. The dosing limit region 2002 includes options for setting a limit on doses per hour, a timed medication delivery limit, or other limits. A user selects one of the options for setting the dosage limit.

The patient specific parameter region 2004 sets parameters related to the dosage limits coordinated to the options in the dosing limit region 2002. The patient specific parameter region 1804 includes a timed delivery limit region 2012, an hourly demand doses region 2014, and a reservoir region 2016.

The timed delivery limit region 2012 sets the delivery limit on a per assigned time period when the option for setting a limit on doses per hour is selected in the dosing limit region 2002. The timed delivery limit region 2012 includes a meter, shown as a slider bar 2018 and indicator 2019. The timed delivery limit region 2012 also includes an initial value gauge 2020, hard limit gauges 2022, soft limit gauges 2024, and control options 2026a-2026c. Operation of the slider bar 2018, indicator 2019, gauges 2020-2024, and control options 2026a-2026c is analogous to operation of the slider bar features as discussed in conjunction with FIG. 18, above.

The hourly demand doses region 2014 sets the delivery limit on a per hour basis when the doses per hour limit is selected in the dosing limit region 2002. The maximum demand doses region 2014 includes a slider bar 2028, indicator 2029, gauges 2030-2034, and control options 2036a-2036c, operation of which is likewise analogous to operation of the slider bar features discussed in FIG. 18.

The timed delivery limit region 2012 and hourly demand doses region 2014 are operated in the alternative, in that only one of the two regions is active at one time. The region that is active depends upon the option selected in the dosing limit region 2002. Selection of a timed delivery limit in the dosing limit
5 region 2002 activates the timed delivery limit region 2012 and deactivates the hourly demand doses region 2014. Selection of an hourly delivery limit activates the hourly demand doses region 2014 and deactivates the timed delivery limit region 2012.

The reservoir region 2016 sets the initial volume settings and display
10 settings for tracking the volume of fluid in the reservoir attached to the medical infusion pump. The reservoir region 2016 includes a meter, shown as a slider bar 2040 and indicator 2041, operation of which is analogous to the slider bar and indicator discussed in conjunction with FIG. 18. The reservoir region further includes an initial value gauge 2042, a disable option 2044, and control options
15 2046a-2046b. The initial value gauge 2042 provides an alternate method for setting the initial value of the reservoir volume to the slider bar 2040 and indicator 2041. The disable option 2044 disables the reservoir volume monitor. The control options 2046a-2046b provide display and edit capabilities to a user of the pump.

The remaining regions, i.e. the titration region 2006, the maximum
20 deliver rate gauge 2008, and the maximum clinician bolus gauge 2010, set pump specific settings related to drug delivery limits. The titration region 2006 enables or disables titration in the medical infusion pump, and sets an optional titration limit in the pump. The maximum delivery rate gauge 2008 sets a maximum delivery rate for the infusion pump. The maximum delivery rate is measured in milliliters per hour,
25 and includes both the basal delivery rate and bolus delivery. The maximum clinician bolus gauge 2010 sets the maximum bolus which can be delivered by a caregiver. The maximum clinician bolus may be a larger bolus than the standard patient-controlled bolus, but must be administered under the supervision of a caregiver. Other regions can be included in the tab 1820 as well.

30 FIG. 21 shows the parameter user interface 1800 with the alarm tab 1830 selected. The alarm tab includes a number of regions used for maintenance and hardware-related alarms, as opposed to the drug delivery threshold alarms discussed above in conjunction with the slider bars. The alarm tab 1830 allows the

user to enable and disable alarms in the medical infusion pump. The alarm tab 1830 includes a reservoir low alarm region 2102, a reservoir empty alarm region 2104, a pump stop alarm 2106, a maintenance alarm 2108, and a hardware alarm 2110. The reservoir low alarm region 2102 provides an alarm indicator when a drug supply
5 volume falls below a threshold level. The threshold level can be a standard level or an administrative user-set volume level. The reservoir empty alarm region 2104 sets a single occurrence or repeating alarm when a drug supply volume falls below a threshold level. The pump stop alarm 2106 sets an alarm which occurs when the medical infusion pump stops operating. The maintenance alarm 2108 enables a
10 maintenance alarm, which alerts a user when maintenance is needed. The hardware alarm 2110 provides options for detection of optional components used with the infusion pump. For example, the hardware alarm 2110 can trigger upon detection of an air detector or other component. The pump hardware detector region 2110 can provide the option of enabling an alarm sent to a caregiver.

15 FIG. 22 shows the parameter user interface 1800 with the security tab 1840 selected. The security tab 1840 includes options related to security of the medical infusion pumps. The security tab 1840 includes a security region 2202 and an epidural region 2204. The security region 2202 presents a number of security options for controlling access to each medical infusion pump. Security options can
20 include an automatic lock, a lock code, clinician code, and an initial lock level. The epidural region 2204 selectably places the pump into a mode configured for epidural therapy.

FIG. 23 shows the parameter user interface 1800 with the display/sound tab 1850 selected. The display/sound tab 1850 includes options
25 related to the display and sound settings for each medical infusion pump. The display/sound tab 1850 includes an auto-review option 2302, a main display region 2304, a units option 2306, a date format option 2308, a send protocol region 2310, and a sound region 2312. The auto-review option 2302 enables a power-up display of options in the medical infusion pump. The main display region 2304 presents a
30 number of programmable fields that will, by default, be displayed on the medical infusion pump. For example, the main display region 2304 includes programmable text display entry, power source information, and drug delivery rate. Other display options related to properties of the medical infusion pump are available as well. The

units option 2306 selectably assigns a units location for display on the medical infusion pump. The date format option 2308 assigns date formats for displaying on the medical infusion pump. The send protocol region 2310 sets one or more patient markers or date/time stamps in the infusion pump network upon distribution of programming instructions to the medical infusion pump. The sound region 2312 enables sound in the medical infusion pump, such as beeping sounds when one or more keys/key sequences are depressed.

FIG. 24 shows the parameter user interface 1800 with the report tab 1860 selected. The report tab 1860 presents options for displaying reports by a medical infusion pump related to events occurring in the pump. The report tab 1860 includes a new patient marker region 2402 and a custom report region 2404. The new patient marker region 2402 includes a list of options to perform related to report generation upon association of a medical infusion pump with a new patient. For example, the new patient marker region 2402 includes a number of report-clearing options and power-on options to be performed when a medical infusion pump is assigned to a new patient. The custom report region 2404 provides a number of selectable options for display in a custom report regarding events tracked in the medical infusion pump network. The custom report region 2404 provides all of the options tracked in the infusion pump network 500, which can be stored in the log files 516. A user selects one or more of the options to generate a report. For example, the custom report region 2404 can include dose counters, doses per hour, pain scale information, drug delivery information, an event log, and a patient marker. The report generated by the custom report region 2404 options displays on the screen of the medical infusion pump 102 or computing system, 104. Both the new patient marker region 2402 and the custom report region 2404 can include additional options as well.

Referring now to FIG. 25, a library export screen 2500 is optionally used to save the pump parameters and protocols for exporting from the database 504 of FIG. 5. One or more libraries are exported using the library export screen 2500, including all of the pump protocols, and containing all of the defined therapies, qualifiers, drugs, and combinations thereof. The library export screen 2500 corresponds to the export library module 724 of FIG. 7, and provides file extraction, in that one or more protocol libraries built using the administrative software 700 are

exported to a data file or database. The export library screen 2500 includes a library field 2502 and an export option region 2504. The library field 2502 displays a number of libraries available to the system that can be exported to a portable file. The export option region 2504 provides export options to the administrative user seeking to export the library data to a file. For example, export options can include creating a read-only file, or a read-write template for use in another instantiation of the software system disclosed herein. The read-only file might be selected in the case that the library is to be loaded onto a medical infusion pump that is disconnected from the database 504 of FIG. 5. The read-write file might be selected if the library is to be transferred to a separate infusion pump network 500 altogether, in which administrative software on that subsequent infusion pump network may be used to subsequently edit that library. Additionally, the export option region 2504 includes an assignable password to add security to the library file exported such that a user attempting to access the protocols contained in the library file must know and enter the correct password.

The above description and figures corresponding to the administrative software 700 provides a therapy-centric programming schema for a medical infusion pump. For example, a certain drug used in conjunction with a medical infusion pump may be appropriate for use with a specific therapy for an adult, but may not be appropriate for the same therapy for a child. Certain drugs may only be appropriate in certain therapies, and under certain qualifying conditions. Pump parameters are initially set according to the protocols defined in the administrative software 700, but are customizable on a pump-by-pump basis using user software associated with a specific pump and/or patient.

FIG. 26 illustrates exemplary architecture of user software 2600 for accessing a pump application program and programming a medical infusion pump. The software 2600 can operate within the pump 102, computing system 104, or a combination thereof. The user software 2600 allows a user, for example a doctor, nurse, pharmacist, or other caregiver, to select and customize pump application protocols, and parameters for execution in and control of a medical infusion pump 102. Depending upon the pump configuration, the user may select a protocol or a library for loading into a medical infusion pump. Additional data structures could be loaded into the medical infusion pump as well. Although the user software 2600

is discussed in conjunction with the administrative software 700 previously described in FIGS. 7-25, it is understood that the user software systems described herein are operable in conjunction with additional hardware/software embodiments.

The medical infusion pumps as described store pump data in
5 memory, such as the memory shown above in FIG. 3. The pump data can include pump parameters, parameter values, programs, and other functional and data systems configured to operate the medical infusion pump. As referred to herein, a set of pump parameters can include the entire memory contents of the pump, or can include a subset of the memory contents, such as selected data values, that can be
10 altered to change operation of the pump.

The user software 2600 is instantiated by a start module 2602. The start module 2602 corresponds to initial execution of the user software 2600 by clicking on an icon on the computer or by some other mechanism for executing software. Upon execution of the start module 2602, the user software 2600 connects
15 to a library on a server 206 containing one or more pump protocols.

Following the start module 2602, operational flow optionally proceeds to a library import module 2604. The library import module 2606 provides the ability to import one or more libraries into the software 2600. This feature can be used by a computing system 104 or medical infusion pump 102 that is not
20 connected to the same medical infusion pump network 500 as the database storing the library 504. If the computing system 104 or medical infusion pump 102 is connected to the medical infusion pump network 500, each component is by default connected to the server 206 and database 504.

The libraries available to be imported include pump protocols and
25 parameters, and may have been created using the administrative software of FIGS. 7-25. The collection of pump protocols can be accessed from the server 206 or in one or more individual computing systems 102. The library import module 2604 allows the user to select one or more pump application programs for downloading to a medical infusion pump 102.

30 Once connected to the desired library either by default or via the library import module 2604, operational flow proceeds to a login module 2606. The login module 2606 regulates user rights in the software 2600 by controlling access to the libraries 508 in the database 504 of FIG. 5. User rights define access levels in

the software for users such as doctors, nurses, other caregivers, or patients. A user will have a set access level allowing the user to view or edit pump application protocols and parameters within the user software 2600. Access levels are set using the user rights module 716 of the administrative software 700, described above in
5 conjunction with FIG. 7. Access levels can be set by a user of the administrative software 700 according to a variety of criteria, such as the type of caregiver (e.g., physician, nurse, or pharmacist).

Different access levels also can provide different rights with respect to pump operational parameters. For example, one access level might give a user a
10 right to edit patient specific pump parameters. One access level might permit a user the right to only view and download the patient specific pump operational parameters. Different embodiments can include the ability to provide an access level for a user any combination of rights to edit, view, and/or download pump operational parameters, protocols, or libraries.

15 Once the user is logged in, the user selectively executes three different modules, a protocol selection module 2608, a task module 2610, and a report module 2612.

The protocol selection module 2608 selects a protocol for use with a medical infusion pump from the protocols loaded in the user software 2600. The
20 protocol selection module 2608 guides the user through selection of a therapy, qualifier, and drug combination defined to be a protocol by the administrative software. The protocol selection module 2608 includes a therapy selection module 2614, a qualifier selection module 2616, and a drug selection module 2618 for this purpose. The therapy selection module 2416 selects a therapy to be administered by
25 the drug infusion pump 102. The therapy is one of the therapies included in the library selected in the library import module 2606. The qualifier selection module 2616 selects a qualifier from those associated with the therapy in the library. The drug selection module 2618 selects a drug associated with the therapy and drug. The protocol selection module 2608 further allows customization of the protocol by
30 allowing a user to modify pump parameters, such as the drug delivery rate, the demand dose, the demand dose lockout, drug delivery limits, and reservoir volume.

The task module 2610 guides a user through maintenance and monitoring tasks that are required for each medical infusion pump 102. These

maintenance and monitoring tasks can include pump settings comparison and testing, as well as changing the reservoir holding the drug delivered by the medical infusion pump 102. The task module 2610 includes a pump settings module 2620, a comparison module 2622, and a reservoir module 2624. The pump settings module 2620 compares the local pump settings to a standing order set for the pump, for example by a caregiver who programmed or customized the medical infusion pump. During operation of the pump settings module 2620, the user software 2600 receives a GUID specific to a protocol and generated by the server 206, and stores the GUID on the pump or computing system. The GUID generated by the server and stored on the pump or computing system is made available to the server when the pump or computing system accesses the library to look up and verify the future access of the correct protocol and/or library. This ensures that the pump settings are compared to the correct protocol stored on the server. The comparison module 2622 compares the local pump settings to a protocol, such as the protocols defined using the administrative software 700. The reservoir module 2624 determines if a drug reservoir is nearly empty and guides a patient or caregiver through the drug cartridge changing process occasionally required during use of a medical infusion pump.

The report module 2612 generates a report from preexisting logged information for a selected medical infusion pump 102. The report module includes a location module 2626 and a type module 2628. The location module 2626 requests the location of the pump from which a report is generated, such as a specific pump or a previously saved report stored on the pump or computing system. The type module 2628 presents a number of types of reports which can be generated from the logged information, such as a drug delivery report or event log, and can display the report responsive to a user request.

Operation of the software terminates at an end module 2630. The end module 2630 corresponds to termination of the administrative software 2600 by clicking on a close window button on the computer or by some other mechanism for terminating execution of software.

FIG. 27 shows a library import screen 2700 used to optionally load a library of pump protocols into the user software 2600. The library import screen 2700 corresponds to the login module 2604, and presents a location field 2702 and a library field 2704 used to browse to and select a library from a library database or

file. An optional password field 2706 accepts user input to allow the user to enter a password for accessing a password protected library file, such as one created using the export library screen shown in FIG. 25. Control buttons 2708a-2708b provide confirmation and cancellation options to the user, allowing the user to complete the library access process.

FIGS. 28-32 illustrate an exemplary process and user interface through which the protocol selection module 2608 leads a user to select a protocol for use with a medical infusion pump. The methods, systems and user interfaces described in conjunction with FIGS. 28-32 can be performed either on a computing system, such as a computing system 104 associated with a medical infusion pump 102, or on a medical fusion pump 102 configured to accept a library of pump protocols directly. Medical infusion pumps 102 accepting a single pump protocol use a corresponding computing system 104.

Referring to FIG. 28, operation of the protocol selection module 2608 is instantiated by a start module 2802. The start module 2802 corresponds to initial execution of the user software 2600, or initial selection of the protocol selection module 2608 by clicking on an icon or tab on the computer or by some other mechanism for executing software.

Following the start module 2802, operational flow proceeds to a load protocols module 2804. The load protocols module 2804 populates the user software 2600 with the protocols from the loaded library. For example, the load protocols module 2804 can populate a listing of protocols for selection using the user software 2600.

Following the load protocols module 2804, operational flow proceeds to a select protocol module 2806. The select protocol module 2806 selects a protocol from among the protocols loaded into the user software 2600 by guiding a user through selection of a therapy, qualifier, and drug defining one of the protocols loaded in the software 2600. The select protocol module 2806 corresponds to the therapy selection module 2614, qualifier selection module 2616, and drug selection module 2618 shown in FIG. 26.

Following the select protocol module 2806, operational flow proceeds to a settings module 2808. The settings module 2808 provides editing and customization of the pump parameters assigned to a medical infusion pump as

dictated by the protocol selected by the user. The parameters include, for example, the drug delivery rate, demand dose rate, or demand dose lockout.

Following the settings module 2808, operational flow proceeds to an optional pump programming module 2810. The pump programming module 2810
5 programs a medical infusion pump with the settings both as defined by the protocol and selected by the select protocol module 2806, and as customized by the settings module 2808. The pump programming module 2810 executes if the user software 2600 resides on a computing system 104 connected to a medical infusion pump 102. The pump programming module may not execute if the software 2600 resides on the
10 medical infusion pump 102 itself, because the protocols are already loaded into the pump alongside the library accessed by the user software 2600.

After the pump programming module 2810 completes, operation of the protocol selection module 2608 terminates at an end module 2812. The end module 2812 corresponds to successful programming of the medical infusion pump.

15 FIGS. 29-32 illustrate a user interface 2900 used to guide a caregiver through the pump programming process. The user interface 2900 operates on a computing system associated to a medical infusion pump. The user interface 2900 includes a login button 2902, a connection status indicator 2904, a library indicator 2906, and a protocol tab 2920, a tasks tab 2940, and a reports tab 2960.

20 The login button 2902, when selected, generates a login screen that checks whether a user has the right to access pump programs, protocols, or parameters. The connection status indicator 2904 displays the connection status of the user software. Connection status can include a connection to a medical infusion pump or connection to a networked server. The library indicator 2906 displays the
25 current library loaded using the library import screen 2700 of FIG. 27.

Referring now to FIG. 29, the user interface 2900 is shown with the protocol tab 2920 selected. FIG. 29 corresponds to an initial state of the protocol selection module 2608 following the load protocols module 2804 of FIG. 28. The protocol tab 2920 guides a user through the process of selecting a protocol,
30 customizing one or more parameters in the protocol, and programming a medical infusion pump with the customized pump program. The protocol tab includes a therapy selection field 2908, a therapy notes field 2910, a qualifier selection field

2912, a qualifier notes field 2914, a drug selection field, 2916, a drug notes field 2918, and a continue button 2922.

The therapy selection field 2908 lists the therapies included in the currently loaded library. For example, the two therapies shown are "Epidural" and "Patient Controlled Analgesia". The therapy notes field displays the notes associated with the selected therapy. In the initial state, the therapy selection field 2908 and therapy notes field 2910 are active, and the qualifier fields 2912, 2914, drug fields 2916, 2918, and continue button 2922 are inactive. No therapy is initially selected in the therapy listing field 2908, so the therapy notes field 2910 remains empty.

FIG. 30 shows the user interface 2900 with the protocol tab 2920 selected and a therapy selected from the therapy selection field 2908. Notes related to the selected therapy appear in the therapy notes field 2910, and the qualifier selection field 2912 and qualifier notes field 2914 activate. A listing of qualifiers associated with the selected therapy appears in the qualifier selection field 2912. The therapy notes shown recite "Epidural Patient Controlled Analgesia" corresponding to the selected therapy, but could contain particular information related to the therapy, such as warnings, descriptions, or other information about application of the therapy. The qualifiers, which appear once the therapy is selected, are shown to include "Adult and Child over 5" and "Child 5 years and under". No qualifier is initially selected, so the qualifier notes field 2914 remains empty. The drug selection field 2916, drug notes field 2918, and the continue button 2922 remain inactive.

FIG. 31 shows the user interface 2900 with the protocol tab 2920 selected and both a therapy selected from the therapy selection field 2908 and a qualifier selected in the qualifier selection field 2912. Notes related to the qualifier appear in the qualifier notes field 2914, and the drug selection field 2916 and drug notes field 2918 are active. For example, "Adult and Child over 5" is shown to be, and the qualifier notes field 2914 displays specific notes applicable to those patients. A listing of drugs associated with the therapy and qualifier appears in the drug selection field 2916. Three exemplary drug menu listings including "Fentanyl 10 mcg/ml", "HYDRO Morphine 1 mg/ml" and "Morphine 1 mg/ml" are shown. No

drug is initially selected, so the drug notes field 2918 remains empty. The continue button 2922 remains inactive.

FIG. 32 shows the user interface 2900 with the protocol tab 2920 selected and a therapy, qualifier, and drug selected in each of the respective selection fields 2908, 2912, and 2916. Once a therapy, qualifier, and drug are selected a specific protocol is designated from among the protocols defined in the administrative software 700. Each notes field 2910, 2914, and 2918 displays information related to the selected therapy, qualifier, and drug, respectively. The continue button 2922 activates, allowing the user to continue with customization of the selected protocol by changing one or more pump parameters related to the therapy, qualifier, and drug.

FIG. 33 illustrates an exemplary process by which the administrative software 2600 customizes one or more parameters of a selected pump protocol, and corresponds to the settings module 2608 of FIG. 26. The settings module 2608 allows customization of one or more pump parameters while monitoring whether the customized value is within an acceptably safe dosage or drug delivery range.

The settings module 2608 is instantiated by a start module 3302. The start module 3302 corresponds to selection of the confirmation button 2922 of FIGS. 29-32.

Following the start module 3302, operational flow proceeds to a display module 3304. The display module 3304 displays the default protocol settings of the protocol loaded onto a medical infusion pump screen or a computing system associated with the pump. The display module 3304 presents a number of meters to a user related to the drug delivery rates and other parameters controlled by the pump. The meters provide user controls for modifying one or more pump parameters.

In various embodiments, the display module 3304 can be configured to display a variety of coloring and image features. In one possible embodiment, the meters are slider bars that can include graphical thresholds set by the administrative software. A cautionary color change of the user interface (i.e. green or gray to yellow or red) can represent a warning to the user that the current setting is outside of the administratively set thresholds.

In another possible embodiment, the overall background color of the user interface is color-coded to correspond to hospital coding procedures, and can represent one or more location-specific warning or status conditions. Additionally, the color coding can be located behind an image displayed on the pump screen, and
5 can be keyed to a location of the medical infusion pump, the drug administered, or a warning condition within the pump.

Following the display module 3304, operational flow proceeds to a custom settings module 3306. The custom settings module 3306 receives the current customized pump settings based on user-customization of one or more pump
10 parameters. The custom settings module 3306 can provide a user customization interface for setting pump parameters to values other than the initial or default values set in the administrative software 700.

Following the custom settings module 3306, operational flow depends upon the implementation of the hard limits and soft limits in the
15 administrative software 700. If the hard limit gauges in FIGS. 18-20 provide a warning but do not dictate an absolute maximum/minimum for the range of programmable values within the user software 2600, operational flow proceeds to a hard limit determination operation 3308. If the hard limit gauges dictate an absolute maximum/minimum for the range of programmable values within the user software
20 2600, the hard limit will likely not be exceeded, so operational flow proceeds directly to the soft limit determination operation 3312.

The optional hard limit determination operation 3308 determines if the pump settings are outside the "hard limits" set in the administrative software 700. If the pump settings exceed the hard limit (i.e. above the maximum or below
25 the minimum value), operational flow branches "yes" to a hard limit indicator module 3310. If the pump settings do not exceed the hard limit, operational flow branches "no" to a soft limit determination operation 3312.

The optional hard limit indicator module 3310 executes in conjunction with the hard limit determination operation 3308, and generates an
30 indicator to a user that the hard limit set in the administrative software is exceeded by the current settings of the medical infusion pump. If the hard limit determination operation 3308 is bypassed or otherwise absent from the user software 2600, the hard limit indicator module 3310 can be absent/bypassed as well. The hard limit

indicator module 3310 creates an alert indicator on the display of the pump or associated computing system, or sends an alert to the server or other computing system to alert a caregiver that an alert condition has been reached by the pump due to exceeding the hard limit. Operational flow proceeds to the display module 3704 to update the display and to allow additional user modification of the pump settings.

The soft limit determination operation 3312 determines if the pump settings are outside the "soft limits" set in the administrative software 700. If the pump settings exceed the soft limit, operational flow branches "yes" to a soft limit indicator module 3314. If the pump settings do not exceed the soft limit, operational flow branches "no" to return to the display module 3304.

The soft limit indicator module 3314 generates an indicator to a user that the soft limit set in the administrative software is exceeded by the current parameter settings. The soft limit indicator module 3314 creates an alert indicator different from the hard limit indicator module 3310 if the hard limit indicator module 3310 exists or executes within the software 2600. For example, the soft limit alert indicator can be a different color, display a different message, or send a different alert to a remote medical care provider.

Following the soft limit indicator module 3314, operational flow proceeds to the display module 3304 to update the display and allow additional user modifications of the pump settings. Upon termination of operation of the medical infusion pump, operational flow terminates at the end module 3316.

Referring now to FIGS. 34-35, an exemplary user interface 3400 for customizing pump parameters is shown. The user interface 3400 corresponds to the settings module 2808 of FIG. 28, and operates generally as described in FIG. 33. The user interface 3400 includes a status indicator 3402, a continuous rate region, 3404, a demand dose region 3406, a demand dose lockout region 3408, a timed limit region 3410, and a reservoir region 3412. The user interface 3400 further includes control buttons 3414a-3414c.

The regions 3404-3412 correspond to the patient specific pump parameters 512a of FIG. 5, and can include one or more of continuous rate, demand dose, demand dose lockout, timed limits, reservoir volume, and other patient-specific parameters. Only those patient specific pump parameters that are associated with the selected index of therapy, qualifier, and drug appear in the user interface

3400, and can be as few as one parameter and can incorporate as many parameters as are programmable within a medical infusion pump.

The user interface 3400 presents a standardized interface to a patient or caregiver using the user software 2600, such as at the point of care of a patient or other location on the infusion pump network. The user interface 2600 corresponds to any of a number of types of user software 2600 and administrative software 700. The user interface 3400 also can be configured to be used with various types of medical infusion pumps 102. Most pumps require programming with both patient specific pump parameters and non-patient specific pump parameters, but vary as to the data structure in which this data is passed to the pump. The user interface 3400 reflects from the user software 2600 the number of regions to create corresponding to the number of patient specific pump parameters, which are generic to various types of pumps. Therefore, the user interface 3400 can be used with any of a number of pumps having various brands, interfaces, data structures, or other variances.

The status indicator 3402 displays the library, therapy, qualifier, and drug which define the protocol loaded by the user software 2600. In the exemplary user interface, the library is an "ICU" library, the therapy selected is "Patient Controlled Analgesia", the qualifier is "Adult and Child over 5", and the drug is "Fentanyl 10 mcg/ml".

The continuous rate region 3404 defines the continuous, or basal, rate of drug delivery in the specific medical infusion pump. The continuous rate region 3404 includes a numerical reading 3416 and a meter, shown as a slider bar 3418 and an indicator 3419. The meter generally has two or more locations, each corresponding to a parameter value that can be programmed in the medical infusion pump. Generally, the positional relationship of the meter indicates the setting of the meter. The numerical reading 3416 reflects the current value for the continuous drug delivery rate.

In the embodiment of the meter shown as the slider bar 3418, the indicator 3419 slides along the slider bar 3418, and the positional relationship between the slider bar 3418 and indicator 3419 dictates the continuous drug delivery rate. Threshold indicators 3420 determine the safe limits within which the continuous rate can be set, and can represent either the hard limits or the soft limits

set for the parameter by the administrative software 700. In the embodiment of the user software 2600 whose absolute threshold levels are limited by the hard limits set in the administrative software 700, the ends of each bar 3418 represent the hard limits and the threshold indicators 3420 represent the soft limits for the continuous rate. The thresholds are tested using the method described in FIG. 33.

The demand dose region 3406 customizes the demand dose, or bolus, delivered by the medical infusion pump. The demand dose region includes a numerical reading 3422 and a meter, shown as a slider bar 3424 and indicator 3425. The slider bar 3424 includes threshold indicators 3426 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3422, slider bar 3424, indicator 3425, and threshold indicators 3426 operate analogously to those in the continuous rate region 3404, but set the bolus level parameter rather than the continuous rate parameter.

The demand dose lockout region 3408 customizes the time period after a bolus is delivered in which no additional bolus can be provided. The demand dose lockout region 3408 includes a numerical reading 3428 and a meter, shown as a slider bar 3430 and indicator 3431. The slider bar 3430 includes threshold indicators 3432 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3428, slider bar 3430, indicator 3431, and threshold indicators 3432 operate analogously to those in the continuous rate region 3404, but set the demand dose lockout period parameter rather than the continuous rate parameter.

The timed limit region 3410 customizes the amount of the selected drug deliverable by a medical infusion pump within a specified timeframe. The timed limit region 3410 also includes a numerical reading 3434 and a meter, shown as a slider bar 3436 and indicator 3437. The slider bar 3436 includes threshold indicators 3438 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3434, slider bar 3436, indicator 3437, and threshold indicators 3438 operate analogously to those in the continuous rate region 3404, but set the timed drug delivery threshold parameter rather than the continuous rate parameter.

The reservoir region 3412 defines the size of the reservoir used in conjunction with the medical infusion pump. The size of the reservoir is relevant to

computing drug delivery volumes for the purpose of setting alarms and other indicators for replacing or refilling the reservoir. The reservoir region 3412, like the other regions, includes a numerical reading 3440 and a meter, shown as a slider bar 3442 and indicator 3443. The slider bar 3442 includes threshold indicators 3444
5 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3440, slider bar 3442, and indicator 3443 operate analogously to those in the continuous rate region 3404, but set the reservoir volume parameter rather than the continuous rate parameter. In the embodiment shown, no threshold indicators are included in the reservoir region. This is because the reservoir region
10 3412 is allowed to use the entire operational range of the reservoir, since no hard limits are set in the reservoir volume region 2016 of FIG. 20 in the administrative software. However, in additional embodiments, one or more threshold indicators can be incorporated into the reservoir region to trigger a warning when the drug reservoir associated with the medical infusion pump contains less than the volume of
15 drugs set by the threshold volume.

The meters in each region 3404-3412 may be adjustable in that each of the patient specific pump parameters are adjustable using a meter. The administrative software 700 enables the adjustment of one or more of the meters by allowing adjustment of those patient specific pump parameters using the options
20 displayed on the user interface 1800 of FIGS. 18-20.

The control buttons 3414a-3414c allow a user to send the currently set parameters to the associated medical infusion pump, cancel the parameter customization, or receive help in the process.

FIG. 35 shows the user interface 3400 wherein the indicator 3437 in
25 the timed limit region 3410 resides along the slider bar 3436 at a location outside the range defined by the threshold indicators 3438. A colored region 3502 appears around the slider bar and indicator, providing a visual warning to the user of the user software 2600 that abnormal or unadvisable pump settings exist. The colored region 3502 can change color (i.e. green or gray to yellow or red) indicating that the value
30 is outside a threshold level.

It is noted that additional screen coloration or textual messages can be used to graphically send messages to a user or programmer of the medical infusion pump or associated computing system. For example, a color code system can be

used to reflect a variety of conditions of the medical infusion pump. For example, a color could represent the current coding of the hospital or other health care facility at which the pump may be located. Additionally, the color code can represent a warning condition, a location at which the pump is used, a drug being administered
5 by the pump, or an alert condition. Of course, the color code could represent additional characteristics of the medical infusion pump as well.

The color code can display on the computing system associated with the medical infusion pump, or can be reflected on a monitor associated to the medical infusion pump itself. Text messages can be sent from the server to be
10 displayed on the monitor of the pump or computing system, such as warnings regarding medication, usage tips for the medical infusion pump, or other medical advice. Additionally, the color code can be placed behind images displayed on the pump which can also represent a region of the hospital, an image of the drug being administered, or other background images. Additionally, the screen coloration
15 described can be represented as a flashing screen, a color changing (cross-fading or otherwise) screen, or various other color patterns.

FIG. 36 shows a pump programming user interface 3600 for guiding a user through the process of sending a pump program to a medical infusion pump. The programming screen can display properties of a pump to be programmed, and
20 can include the settings 3602a-3602e that are to be sent to the pump, a protocol indicator field 3604, and a confirmation button 3606.

The settings 3602a-3602e reflect the customized pump parameters set using the user interface 3600 of FIGS. 34-35. The customized pump parameters correspond to the patient specific pump parameters 512a of FIG. 5, and can include
25 one or more of continuous rate, demand dose, demand dose lockout, timed limits, reservoir volume, and other patient-specific parameters.

The protocol indicator field 3604 displays the current protocol selected, in this case shown as "Patient Controlled Analgesia", "Adult and Child over 5", and "Fentanyl 10 mcg/ml". The confirmation button 3606 sends the pump
30 program, including the settings 3602a-3602e for the pump parameters, to the pump.

Once the pump is programmed, the pump program executes according to the protocol selected in conjunction with the customized parameters. Referring now to FIG. 37, an exemplary software process 3700 for displaying

medical infusion pump customizations is shown. The process occurs within the pump 102, computing system 104, or a combination thereof, and can be part of a pump program sent to a medical infusion pump.

Customizations in protocol programming refer to differences between
5 the actual pump operation and a standing order (i.e. settings programmed by an administrative user). The standing order can be an original pump parameter or initial value, and the pump operation for comparison can be either a current pump parameter or simply a non-original pump parameter.

A start module 3702 initiates the process 3700. Operational flow
10 proceeds to a comparison module, shown as a compare pump settings to standing order module 3704. The pump settings stored on the pump or computing system as shown above in FIG. 5 can be compared against a standing order stored on a server as a pump protocol. To ensure that the correct standing order is accessed for comparison, a GUID assigned to the loaded pump settings corresponds to the
15 protocol stored on the server. A display change bar module 3706 presents an indicator on either the medical infusion pump or associated computing system. In the display change bar module, the original pump parameter, or portion of a standing order, can be juxtaposed against the non-original pump parameter, such as in a table format. One or both of the original and non-original pump parameters can be a
20 standard or customized pump parameter. Pump parameters that are displayed in the display change bar module 3706 include drug delivery rate, drug capacity, remaining capacity of the medical infusion pump, bolus levels or occurrences, alarm occurrences, threshold, and/or frequency, drug delivery periods, or other parameters.

In one possible embodiment, a legend can indicate the meaning of the
25 pump parameters being displayed, and a time/date stamp can display the time at which the original and/or non-original pump parameter was measured. In additional embodiments, differences between the original and non-original pump parameters can be highlighted when displayed, such as by using a color change or other indicator.

30 Operational flow terminates at an end module 3708.

FIG. 38 is an exemplary schematic illustration of a change bar 3800 displayed on a medical infusion pump 102. The change bar 3800 as shown is displayed on the screen 406 of the pump described above in FIG. 4. Alternately, the

change bar 3800 can be displayed on a computing system 104 associated with the medical infusion pump 102.

The change bar 3800 includes a plurality of change bar entries 3802. The change bar entries correspond to pump parameters, and in the figure shown the
5 entries 3802a-3802b correspond to patient specific pump parameters. The change bar can display non patient specific pump parameters as well.

The change bar 3800 can compare any of a number of original and non-original pump parameters. In one embodiment, the change bar 3800 compares the current operation of the pump to the originally programmed operation of the
10 pump. In another embodiment, the change bar 3800 compares the operation of the pump as initially programmed to the suggested programming of the pump based on the original pump protocol. In a further embodiment, the change bar 3800 compares historical activity of the pump to the current pump protocol.

In one embodiment of the change bar 3800, the change bar entries
15 3802a-3802b change color when the difference between the original and the non-original pump parameters exceeds a threshold amount. The threshold amount can be, for example, the soft limits set in the administrative software 700. In a further embodiment, the text can change color when a difference greater than the threshold is detected. The change bar can incorporate additional graphics and images on the
20 display.

Referring now to FIG. 39, user software 2600 is again described in which the user interface 2900 of FIGS. 29-32 is shown with the tasks tab 2940 selected. The tasks tab 2940 includes a tasks region 3902 for selection of one or more task options, of which a comparison is displayed when the user selects the
25 "continue" option 3904 shown. Tasks include operations related to maintenance of the pump once in operation. The tasks region 3902 includes a number of pump comparison options, such as a compare pump settings option 3906, a compare pump settings to protocol option 3908, and a change reservoir option 3910. The compare pump settings option 3906 compares the pump settings to the original protocol from
30 which the pump parameters were based. The compare pump settings to protocol option 3908 compares all pump settings for the protocol selected.

The user software 2600 accesses the protocol loaded on the server 206 to compare the current pump settings to the original or current protocol using

the options 3904, 3906. To accomplish this, it is necessary for the user software 2600 to clarify to the server 206 which protocol is being compared within the database 504 of FIG. 5. The user software 2600, in conjunction with the server 206 uses the globally unique identifier (GUID) described above in FIG. 5 to provide the
5 identifier for corresponding the protocol on the pump 102 to the protocol as stored in the server 206. The GUID can be generated by the server 206 and transmitted alongside the protocol and/or library when transmitted to the computing system 104 or infusion pump 102, as described above.

The change reservoir option 3910 guides a user of the software 2600
10 through changing a drug reservoir used in conjunction with the medical infusion pump.

FIG. 40 shows the user interface 2900 with the reports tab 2960 selected. During or after pump operation, reports including drug delivery or event logs are used to detect the condition of a patient or of the medical infusion pump.
15 The events which can be tracked using reports tab 2960 are those which are available due to being automatically tracked by the medical infusion pump. The reports tab 2960 presents a number of selectable options for generating reports of pump activity, including time/date information, drug delivery information, and other event information. The reports tab 2960 includes source fields 4002 and option
20 regions 4004. The source fields 4002 present a variety of sources from which reports can be drawn. The source fields 4002 can include the medical infusion pump and stored reports saved within the network. The option regions 4004 present a number of options related to the selected source. For example, a report generated directly from a medical infusion pump can be produced based on prescription
25 settings, and event log, a patient history, drug delivery, or a reported pain scale. A report generated from a saved report can be produced by indicating the patient identification, the pump identification, or the report type and date. A view field 4006, upon selection by a user, generates the report based on the source and options selected.

30 FIG. 41 shows a report user interface 4100 for displaying operation of a medical infusion pump. The report user interface 4100 shows the report generated using the options selected in the report tab 2960. The report shown in the report screen is a drug delivery report, and can be printed, saved, or discarded by the

user. The drug delivery report can include the volume of the drug delivered, as well as the timing of delivery of the drug. Additional attributes of the medical infusion pump can be reported in the drug delivery report as well.

Aspects of the invention described as being carried out by a
5 computing system or otherwise described as a method of control or manipulation of data may be implemented in one or a combination of hardware, firmware, and software. Embodiments of the invention may also be implemented as instructions stored on a machine-readable medium, which may be read and executed by at least one processor to perform the operations described herein. A machine-readable
10 medium may include any mechanism for storing or transmitting information in a form readable by a machine (e.g., a computer). For example, a machine-readable medium may include read-only memory (ROM), random-access memory (RAM), magnetic disc storage media, optical storage media, flash-memory devices, electrical, optical, acoustical or other form of propagated signals (e.g., carrier waves,
15 infrared signals, digital signals, etc.), and others.

In the foregoing detailed description, various features are occasionally grouped together in a single embodiment for the purpose of streamlining the disclosure. This method of disclosure is not to be interpreted as reflecting an intention that the claimed embodiments of the subject matter require
20 more features than are expressly recited in each claim. Rather, as the following claims reflect, inventive subject matter lies in less than all features of a single disclosed embodiment. Thus, the following claims are hereby incorporated into the detailed description, with each claim standing on its own as a separate preferred embodiment. Therefore, the spirit and scope of the appended claims should not be
25 limited to the description of the preferred versions contained herein.

The Claimed Invention Is:

1. An apparatus for indicating a current condition of a medical infusion pump, the apparatus comprising:
 - a memory configured to store two or more color codes, each color
 - 5 code corresponding to a condition of the medical infusion pump and a color;
 - a monitor;
 - a programmable circuit in electrical communication with the memory and the monitor, the programmable circuit programmed to generate a graphical user interface displaying a color upon detection of a condition.
- 10 2. The apparatus of claim 1, wherein the programmable circuit is configured to display the color on a portion of the monitor.
3. The apparatus of claim 1, wherein the condition is selected from the group consisting of:
 - a portion of a hospital coding system;
 - 15 a medical infusion pump warning condition;
 - a location at which the pump is used;
 - a drug currently administered by the medical infusion pump; and
 - an alert condition.
4. The apparatus of claim 1, wherein the color code includes an image
- 20 representing a graphical representation of a region of a hospital.
5. The apparatus of claim 1, wherein the color code includes a background image on the display.
6. The apparatus of claim 1, wherein the programmable circuit is programmed to change the color on the monitor while detecting the condition.
- 25 7. The apparatus of claim 1, wherein the programmable circuit is programmed to cross-fade the color on the monitor while detecting the condition.

8. The apparatus of claim 1, wherein the programmable circuit is programmed to flash the color on the monitor while detecting the condition.
9. The apparatus of claim 1, wherein the color includes patterns of color.
- 5 10 The apparatus of claim 1, wherein the color is a text color.
11. A method of indicating a current condition of a medical infusion pump, the method comprising:
upon detecting a condition in the medical infusion pump, displaying a color on a monitor associated to the medical infusion pump.
- 10 12. The method of claim 11, wherein displaying comprises displaying the color on a portion of the monitor.
13. The method of claim 11, wherein the condition is selected from the group consisting of:
a portion of a hospital coding system;
15 a medical infusion pump warning condition;
a location at which the pump is used;
a drug currently administered by the medical infusion pump; and
an alert condition.
14. The method of claim 11, wherein the color code includes an image
20 representing a graphical representation of a region of a hospital.
15. The method of claim 11, wherein the color code includes a background image on the display.
16. The method of claim 11, further comprising changing the color on the monitor while detecting the condition.

17. The method of claim 11, further comprising cross-fading the color on the monitor while detecting the condition.

18. The method of claim 11, further comprising flashing the color on the monitor while detecting the condition.

5 19. The method of claim 11, wherein the color includes patterns of color.

20. The method of claim 11, wherein the color is a text color.

21. The method of claim 11, wherein the color relates to a difference between an original pump parameter and a non-original pump parameter.

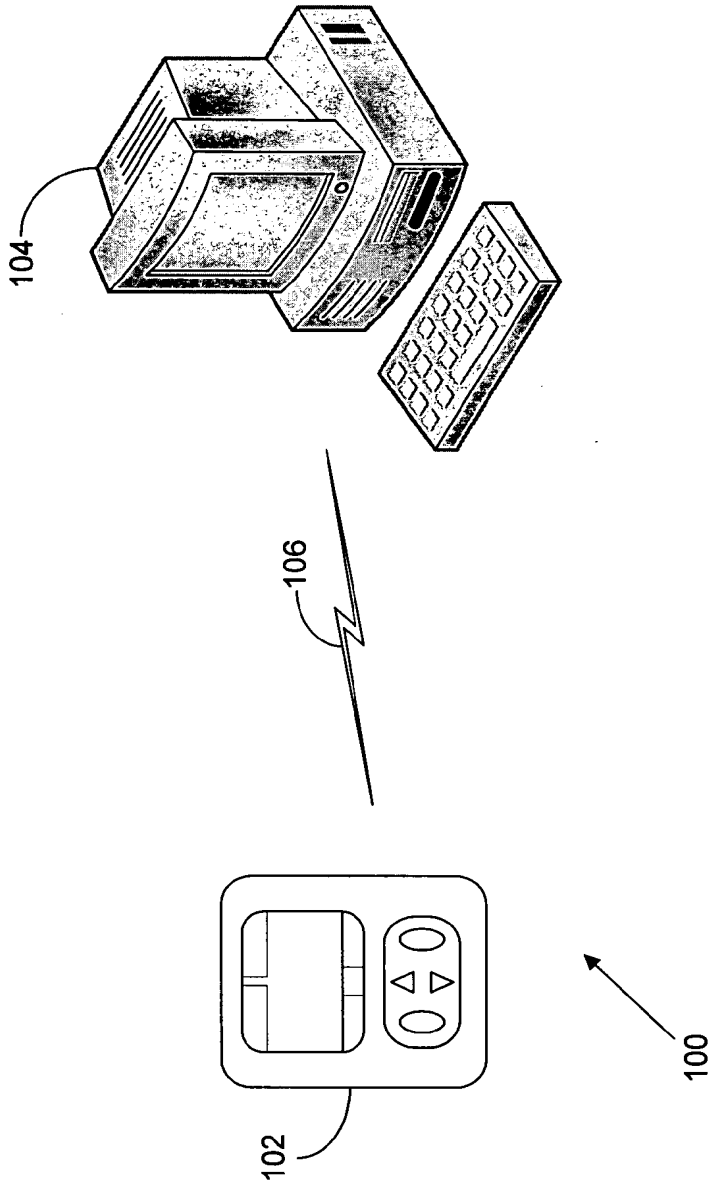


FIG. 1

FIG. 2

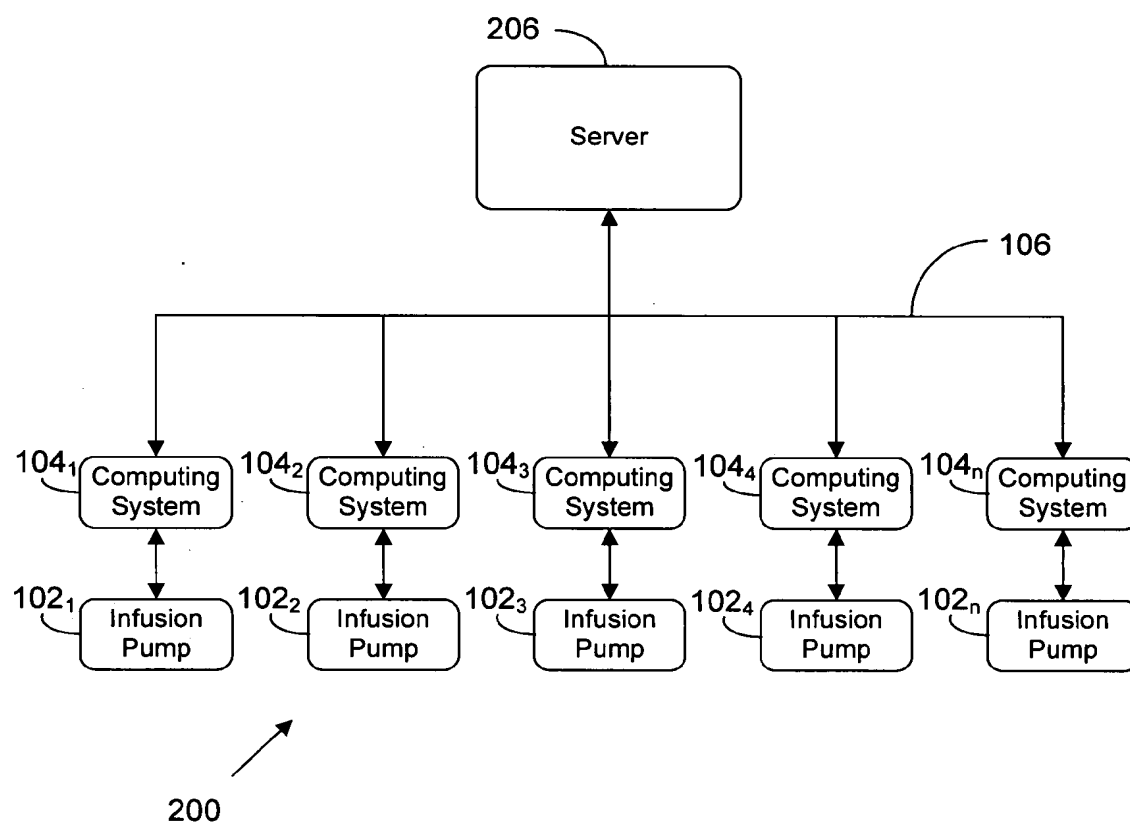


FIG. 3

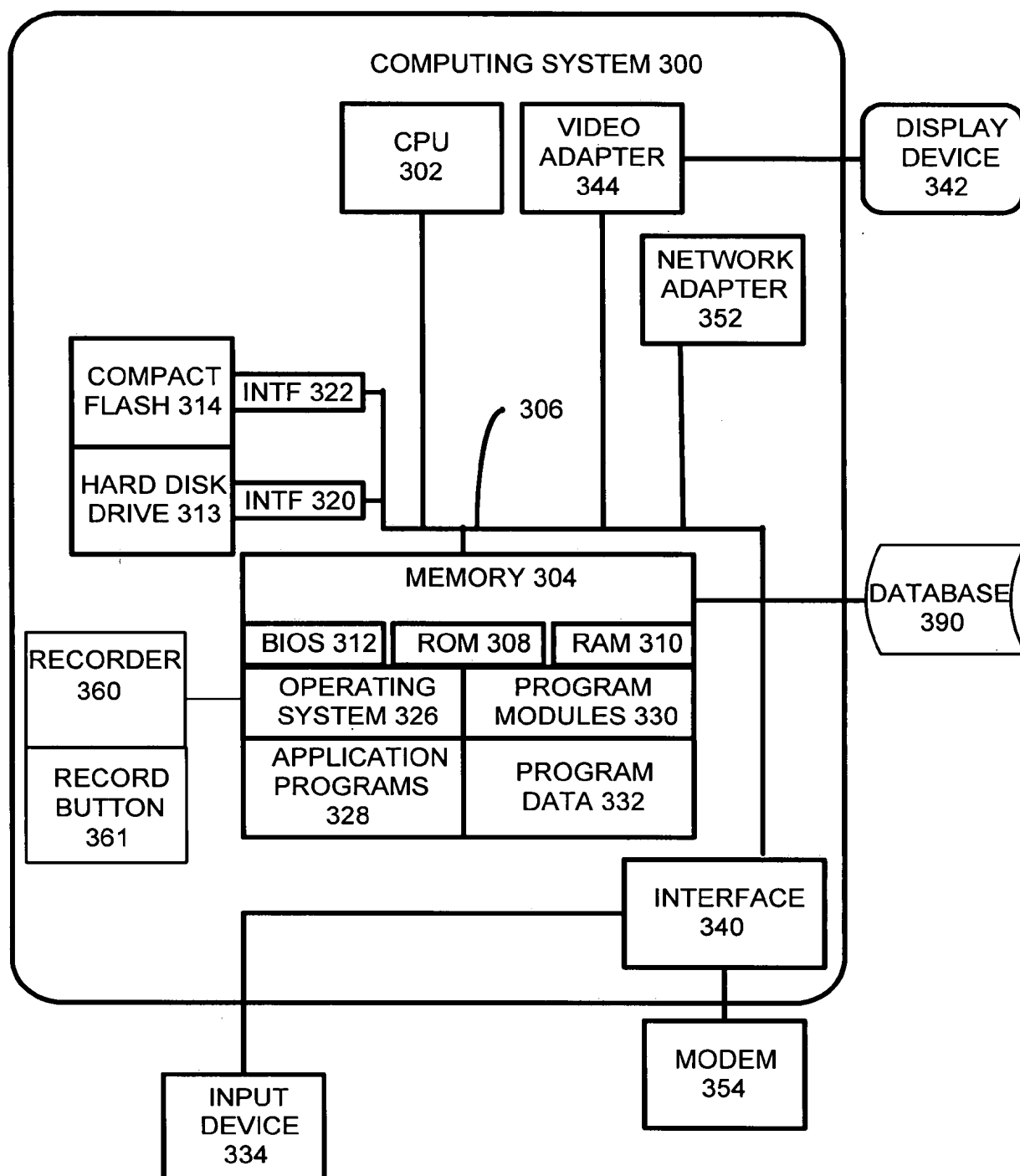


FIG. 4

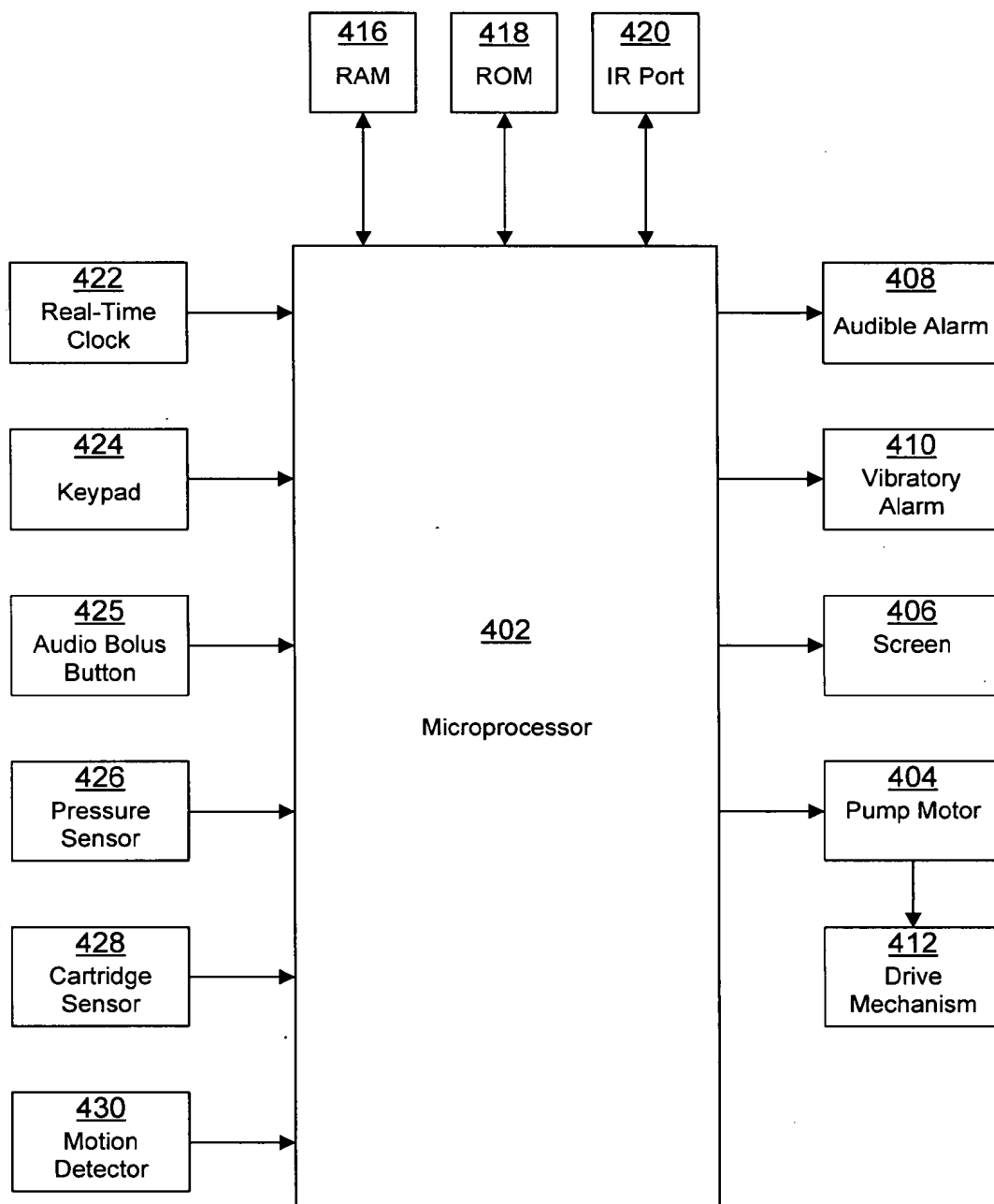
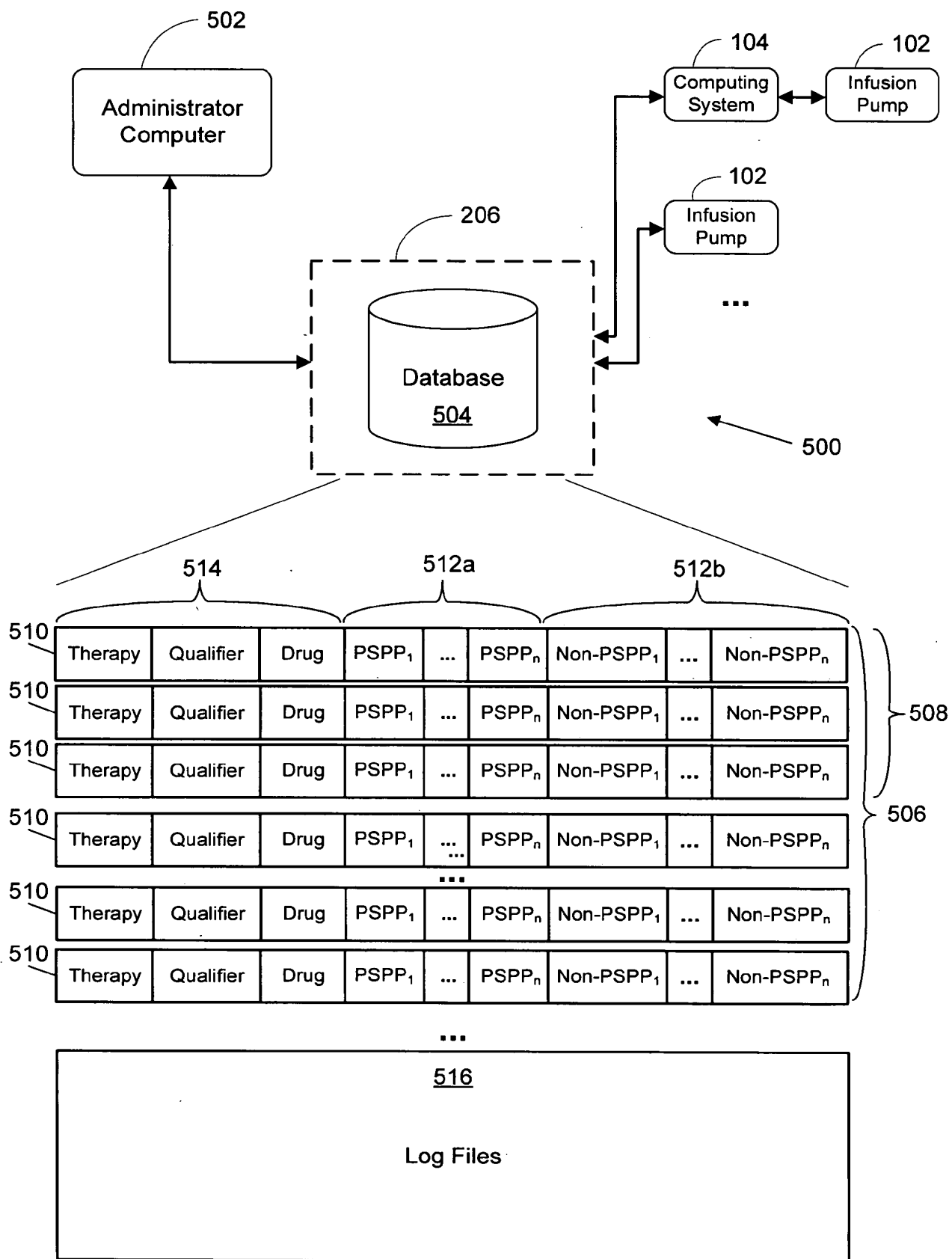


FIG. 5

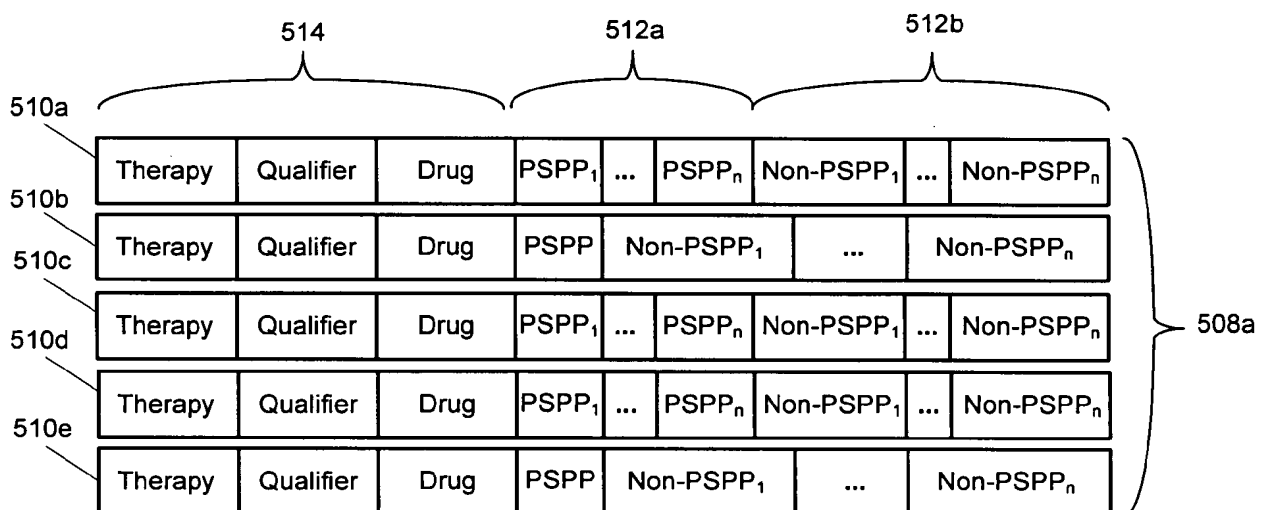


FIG. 6A

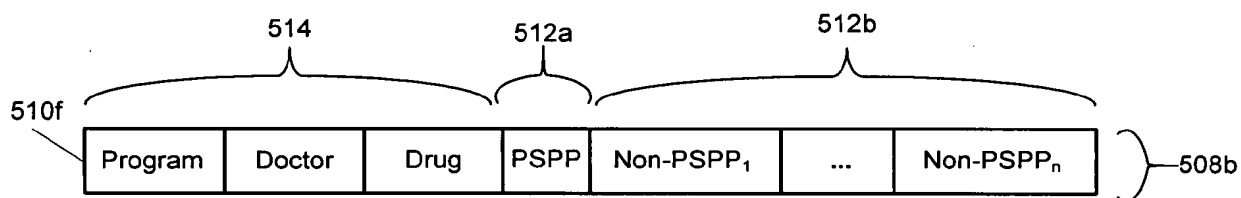


FIG. 6B

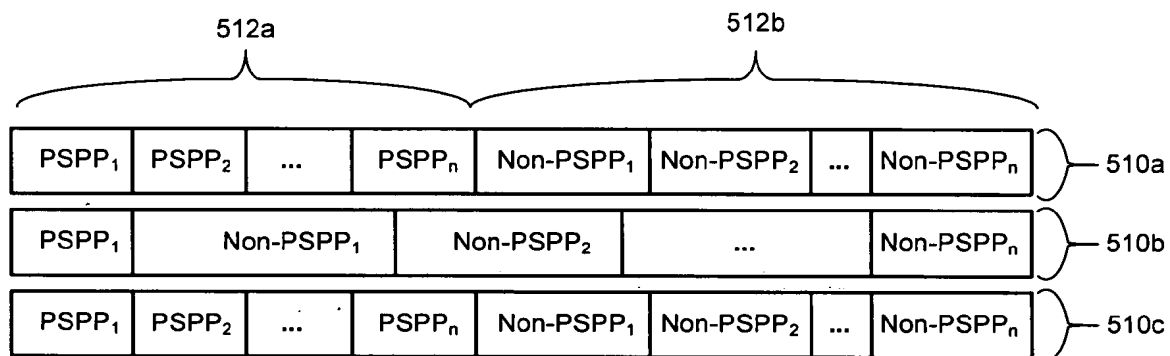


FIG. 6C

FIG. 7

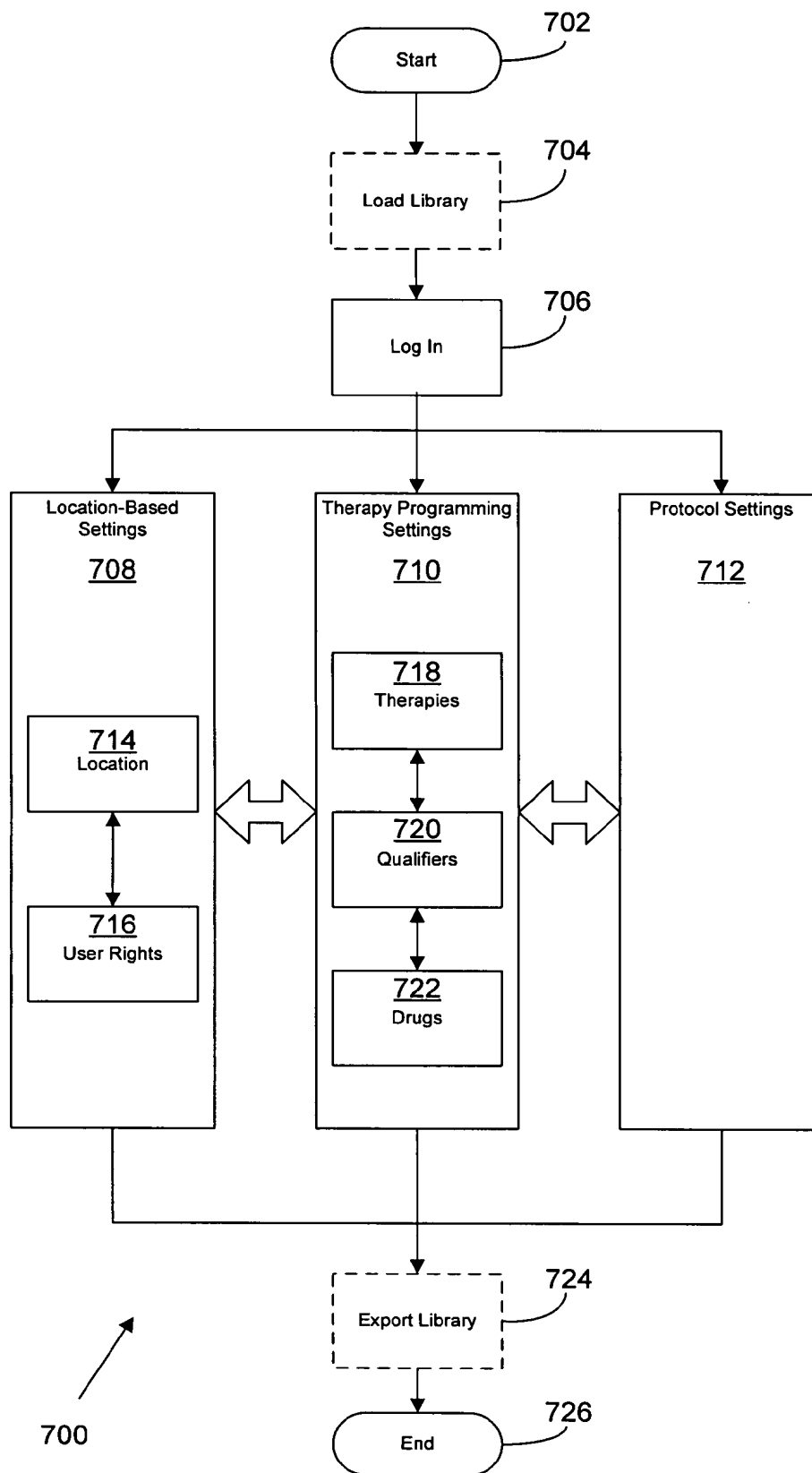


FIG. 8

800

Import Library

Import from Location: 808

Browse 810

Select Protocol Libraries:

- ☐ Neonatal Intensive Care Unit (NICU)
- ☐ Post Anesthesia Care Unit (PACU)
- ☐ Intensive Care Unit (ICU)
- ☐ Our Library Base Template 802
- ☐ Library One
- ☐ Library Two

Select All 806a Clear All 806b

Enter Import Database Password:

804

OK Cancel Help

FIG. 9

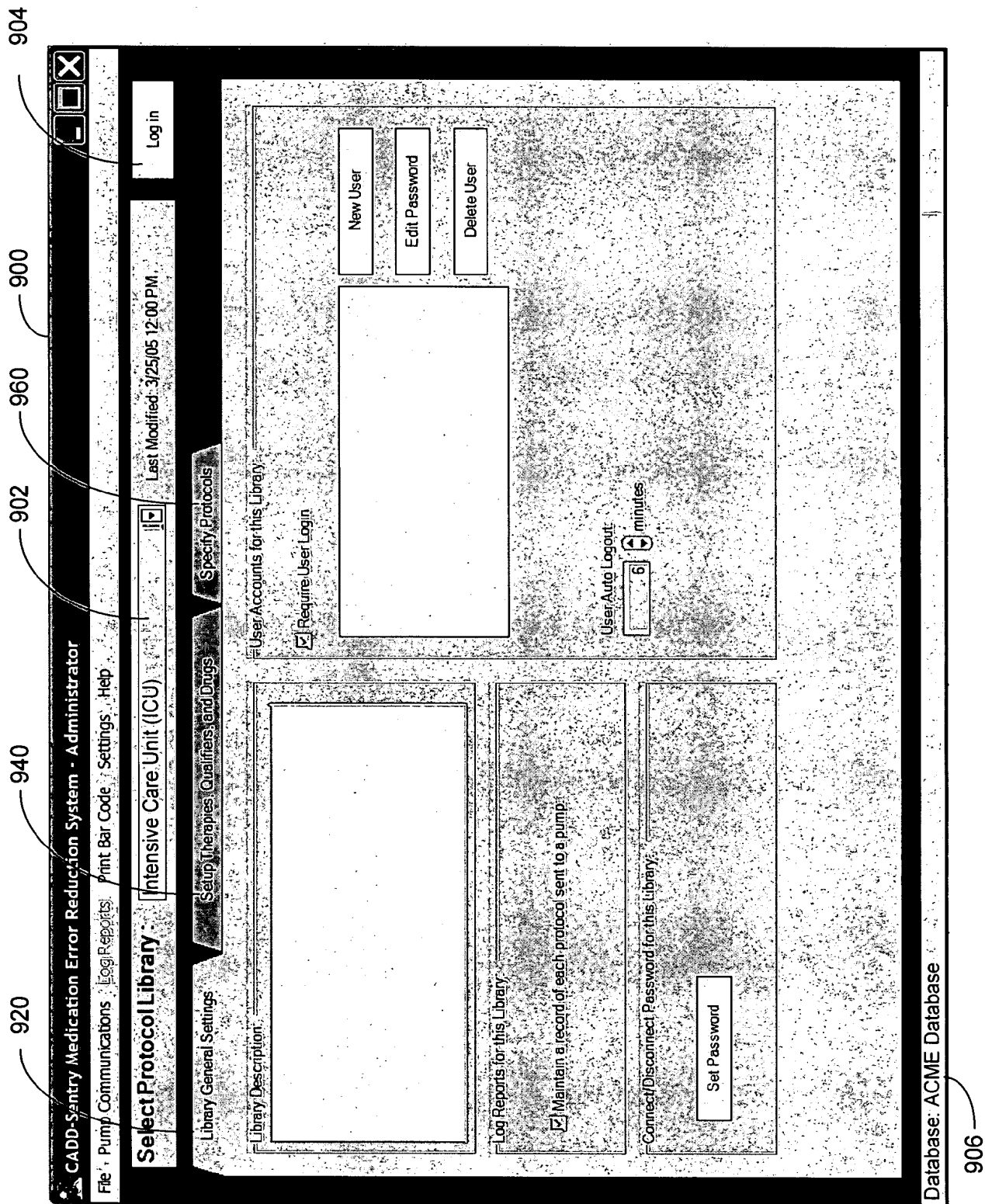


FIG. 10

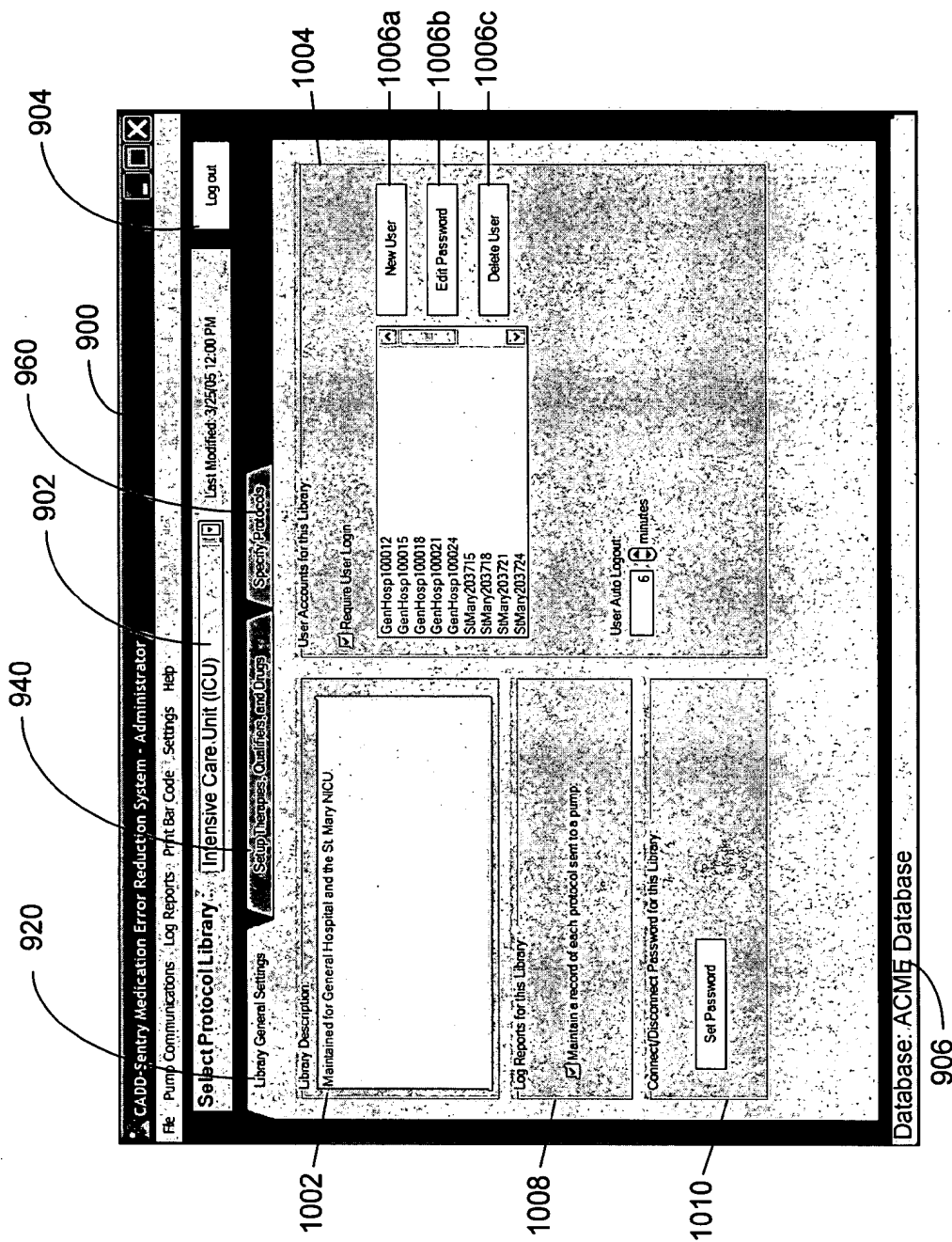


FIG. 11

920 940 960 900 904 1104 1102 1106 1122a 1122b 1122c 1122d 906

CADD-Sentry Medication Error Reduction System - Administrator
 File Pump Communications Log Reports Print Bar Code Settings Help
 Log out

Select ProtocolLibrary: Intensive Care Unit (ICU) Last Modified: 3/25/05 12:00 PM

Library General Settings Setup Therapies Qualifiers and Drugs Specify Protocols

Therapies

Therapy Name	In Use
Patient Controlled Analgesia	Y
Epidural	Y

1108

Notes:
 Patient Controlled Analgesia notes could contain information pertinent to this therapy.

1110 1112a 1112b 1112c New Edit Delete

Qualifiers

Qualifier Name	In Use
Adults over 65 yrs	N
Adults 40 to 65 yrs	Y
Adolescents 10-14 yrs	Y
Children 5-10 yrs	N

1114

Notes:
 Adults over 65 yrs - Additional notes would follow

1116 1118a 1118b 1118c New Edit Delete

Drugs in Neonatal Intensive Care (NICU) Library

Name	Units	Concentration	Drug ID	In Use
Fentanyl	Micrograms	10	123456789	Y
Fentanyl-Morphine	Micrograms	0	esbsssa	N
Hydro-Morphine	Micrograms	30	9999888	N
Hydro-Morphine	Micrograms	80	HydroMorpho	Y
Morphine	Micrograms	0	123456	N

1120

New Drug Edit Drug ID Delete Drug Print Bar Code

Database: AQME Database

FIG. 12

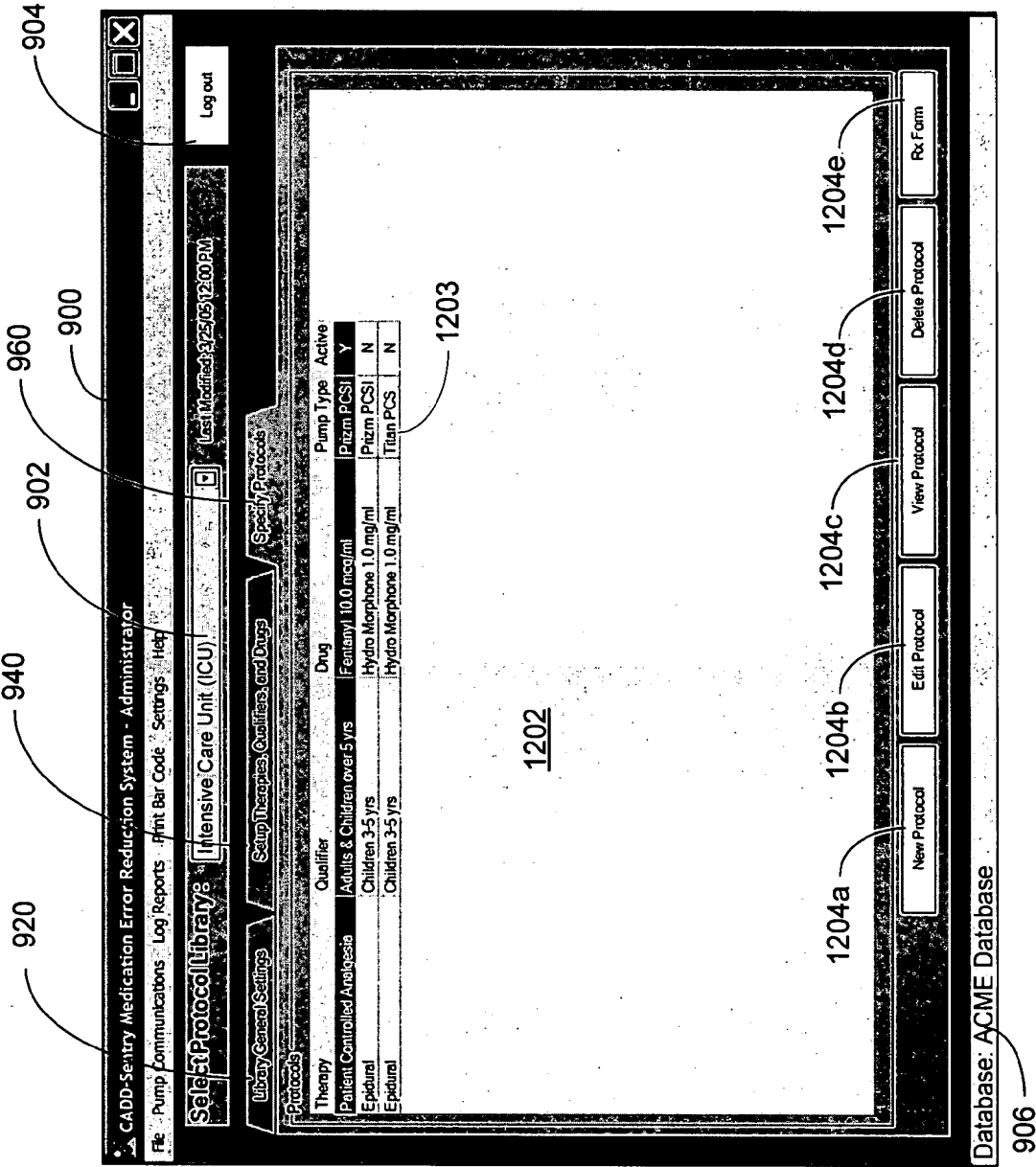


FIG. 13

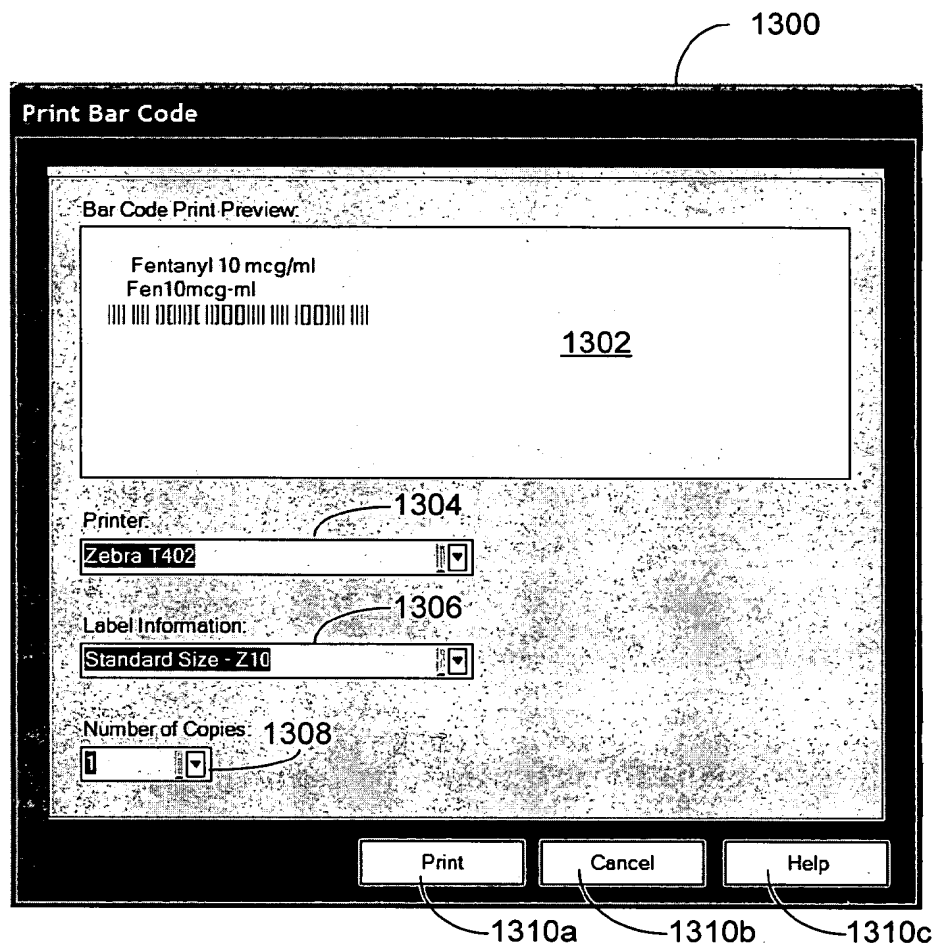


FIG. 14

SUBSTITUTE SHEET (RULE 26)

View Prescription Order Form

77%

1 of 4

<

>

Print

Close

Prescription Form:

Patient Controlled Analgesia
Adults & Children over 5 yrs

Date: _____

Time: _____

Therapy:	Patient-Controlled Analgesia
Notes:	Patient controlled analgesic setup with neo-342 and 688A no 31245.
Qualifier:	Adults & Children over 5 yrs
Notes:	Monitor and record respiratory rate and level of consciousness every two (2) hours times twelve (12) hours upon PCA initiation or dose change, then every four (4) hours. Record pain location and intensity every four (4) hours

Epidural 100.0 mcg/ml

Continuous Rate (mg/hr)	100 to 50 mcg/hr
Demand Dose (mg)	5 to 20 mcg
Dose Lockout (min)	(0.500 to 3:15 hr)
4 Hour Limit (mcg)	(800 to 1000 mcg)
Reservoir Volume (ml)	----
Predicted Maximum Infusion for Side Effect (LSD, etc.) (mcg/hr)	decreased urine output, decreased BP, and vomiting.

Hydro Morphine 1.0 mcg/ml

Continuous Rate (mg/hr)	100 to 50 mcg/hr
Demand Dose (mg)	5 to 20 mcg
Dose Lockout (min)	(0.500 to 3:15 hr)
4 Hour Limit (mcg)	(800 to 1000 mcg)
Reservoir Volume (ml)	----

Predicted Notes: Consider reducing dose gradually. Monitor for side effects such as: decrease anxiety, confusion, decreased urine output, hallucinations, constipation, dry mouth, nausea and vomiting.

Morphine 1.0 mcg/ml

Continuous Rate (mg/hr)	100 to 50 mcg/hr
Demand Dose (mg)	5 to 20 mcg
Dose Lockout (min)	(0.500 to 3:15 hr)
4 Hour Limit (mcg)	(800 to 1000 mcg)
Reservoir Volume (ml)	----

FIG. 15

1500

New Protocol

Therapy Patient Controlled Analgesia	1502	Qualifier 1506	Drug 1510
Patient Controlled Analgesia Epidural			
Therapy Notes Patient controlled analgesia.	1504	Qualifier Notes 1508	Protocol Notes 1512

Next Cancel Help

FIG. 16

1500

New Protocol

Therapy Patient Controlled Analgesia	1502	Qualifier Adults & Children over 5 yrs Children 3-5 yrs	1506	Drug	1510
Therapy Notes Patient controlled analgesia	1504	Qualifier Notes	1508	Protocol Notes	1512

Next Cancel Help

FIG. 17

1500

New Protocol

<p>Therapy</p> <p>1502</p> <p>Patient Controlled Analgesia</p> <p>Therapy Notes:</p> <p>Patient controlled analgesia</p> <p>1504</p>	<p>Qualifier</p> <p>1506</p> <p>Adults & Children over 5 yrs</p> <p>Qualifier Notes</p> <p>Monitor and record respiratory rate and level of consciousness every two (2) hours times twelve (12) hours upon PCA initiation or dose change, then every four (4) hours. Record pain location and intensity every four (4) hours</p> <p>1508</p>	<p>Drug</p> <p>1510</p> <p>Fentanyl 10.0 mcg/ml Hydro Morphine 1.0 mg/ml Morphine 1.0 mg/ml Morphine 10.0 mg/ml</p> <p>1512</p>
--	--	---

Next Cancel Help

FIG. 18

1800

1802

1860

1850

1840

1830

1820

1810

1804

Protocol Details

☐ Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 10.0 mcg/ml

Delivery 1

Delivery 2

Alarms

Security

Display/Sound

Reports

General Settings

☐ Require Drug Verification

1816

☐ Require Second Nurse Verification

(User Logon must be enabled to use this feature.)

Weight-Based

☐ Set as a Weight-Based Protocol

1818

Weight Limit: kg

Limits for Patient Specific Parameters:

Continuous Rate mcg/hr

☐ Display 1822

☐ Allow Edit 1824

☐ Require Comment 1826

Hard Maximum: 1834

Soft Maximum: 1836

Initial Value: 1832

Soft Minimum: 1836

Hard Minimum: 1834

Demand Dose mcg

☐ Display 1844

☐ Allow Edit 1846

☐ Require Comment 1848

Hard Maximum: 1846

Soft Maximum: 1848

Initial Value: 1844

Soft Minimum: 1848

Hard Minimum: 1846

Demand Dose Lockout

☐ Display 1854

☐ Allow Edit 1856

☐ Require Comment 1858

Hard Maximum: 1858

Soft Maximum: 1862

Initial Value: 1856

Soft Minimum: 1862

Hard Minimum: 1858

Save 1806a

Cancel 1806b

Help 1806c

FIG. 19

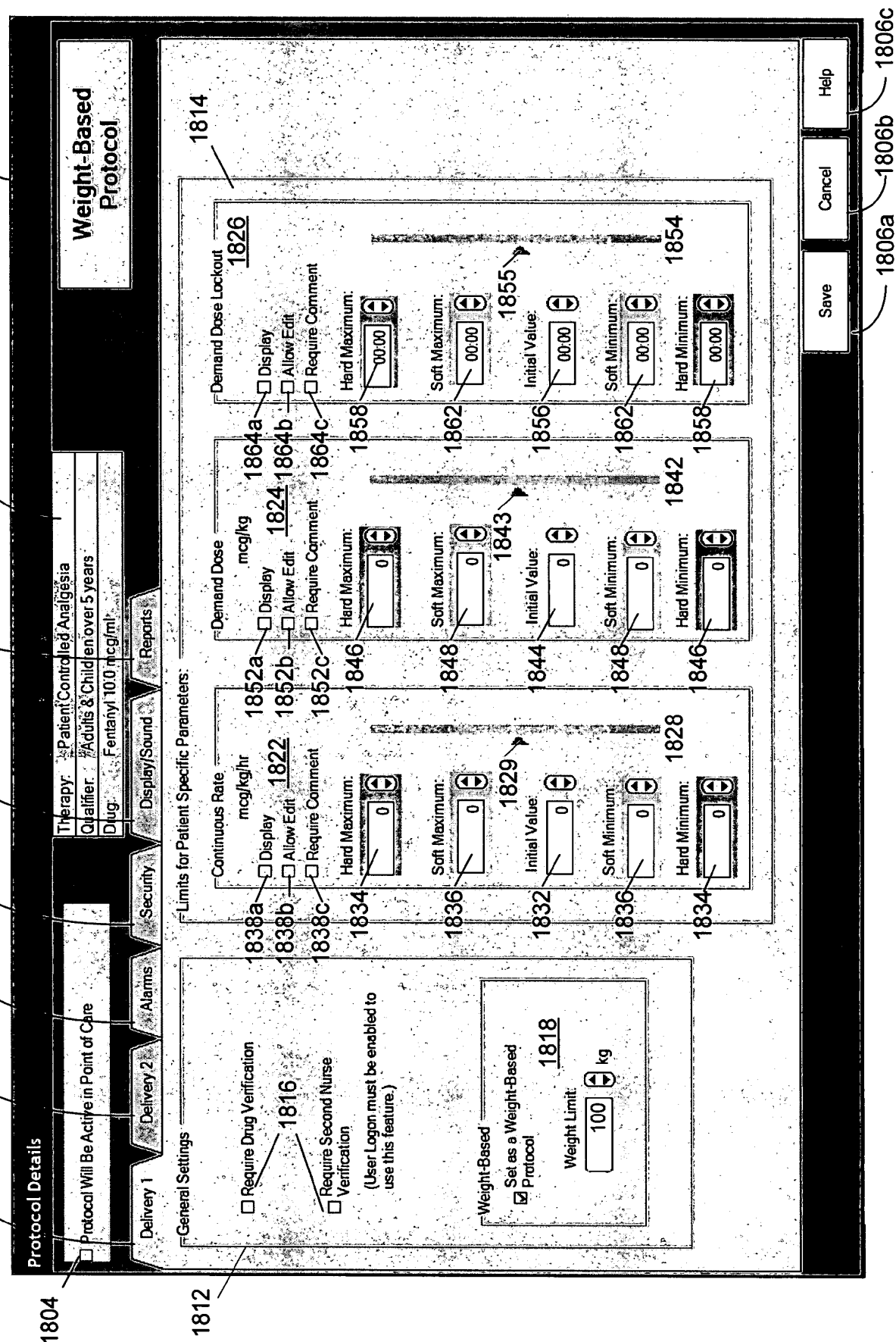


FIG. 20

1800

1802

1860

1850

1840

1830

1820

1810

1804

Protocol Details

Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drip: Fentanyl 10.0 mcg/ml

Delivery 1

Delivery 2

Alarms

Security

Display/Sound

Reports

Limits for Patient Specific Parameters:

Dosing Limit:

☒ Delivery Limit:

☐ No Limit

hours

☐ Doses per Hour

☐ Neither

2026a

2026b

2026c

2022

2024

2020

2024

2022

4 Hour Delivery Limit mcg

2036a

2036b

2036c

2032

2034

2030

2034

2032

2018

2032

2012

2014

2016

Maximum Demand Doses per Hour

☐ Display

☐ Allow Edit

☐ Require Comment

Hard Maximum:

Soft Maximum:

Initial Value:

Soft Minimum:

Hard Minimum:

2042

2044

2046a

2046b

Reservoir Volume ml

☐ Not in Use

☐ Initial Value:

2040

2008

2010

Maximum Delivery Rate: 125 ml/hr

Maximum Clinician Bolus: 0 mcg

2006

Titration

☐ No Titration

☐ Titration Limit

☐ No Limit

2002

2004

2006a

2006b

2006c

Save

Cancel

Help

FIG. 21

1804

1810

1820

1830

1840

1850

1860

1802

1800

Protocol Details

☐ Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 10.0 mcg/ml

Delivery 1

Delivery 2

Alarms

Security

Display/Sound

Reports

Reservoir Low Alarm 2102

☐ Standard

☐ Trip Point: 0 ml

Reservoir Empty Alarm Type 2104

☐ Single

☐ Insistent

Pump Stopped Alarm 2106

☐ Beeps

☐ Two Tone

Preventive Maintenance 2108

☐ Reminder On

☐ Reset Reminder

Interval: 0 months

Due: 03/02/2005

Pump Hardware 2110

☐ Air Detector Required

☐ Air Detector On

☐ Upstream Sensor On

Save 1806a

Cancel 1806b

Help 1806c

FIG. 22

1800

1802

1860

1850

1840

1830

1820

1810

1804

Protocol Details

☐ Protocol Will Be Active in Point of Care

Delivery 1

Delivery 2

Alarms

Security

Display/Sound

Reports

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 100 mcg/ml

Security

Auto Lock:
Not in Use

Lock Code:
61

Clinician Code:
997

Initial Lock Level:
2

2202

Epidural Mode

☐ Epidural Mode On

2204

Save

Cancel

Help

1806a

1806b

1806c

FIG. 23

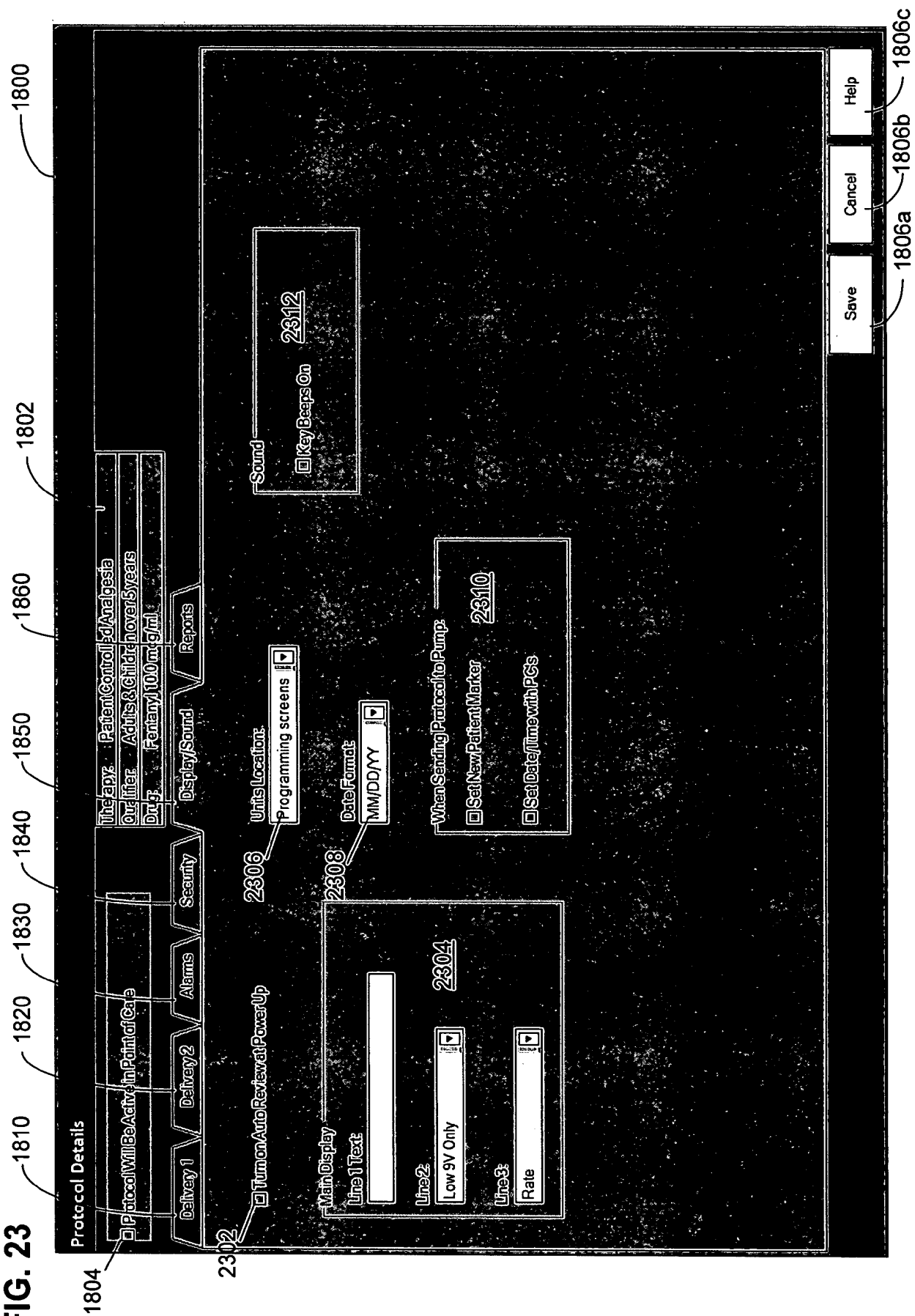


FIG. 24

1804

1810

1820

1830

1840

1850

1860

1802

1800

Protocol Details

☐ Protocol Will Be Active in Point of Care

Delivery 1

Delivery 2

Alarms

Security

Display/Sound

Reports

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 Years

Drug: Fentanyl 10.0 mcg/ml

New Patient Marker

☐ Reports/No Clear

☐ Power Up/No Clear

☐ Reports/Clear

☐ Power Up/Clear

2402

Custom Reports

☐ Dose Counters

☐ Given

☐ Doses Hour by Hour

☐ Patient Review

☐ Pain Scale Entry

☐ Pain Scale Log

☐ Delivery Log

☐ Event Log

☐ New Patient Marker

2404

Save

Cancel

Help

1806a

1806b

1806c

FIG. 25

2500

Export Library

Select Protocol Libraries:

- ☐ Neonatal Intensive Care Unit (NICU)
- ☐ Post Anesthesia Care Unit (PACU)
- ☐ Intensive Care Unit (ICU)
- ☐ Our Library Base Template
- ☐ Library One
- ☐ Library Two

2502

Select All Clear All

Export Option

☒ Read Only for Use at Point of Care

2504

☐ Read-Write Template for Use in Administrator

Create New Passwords:

(This password will be used when importing the selected libraries.)

OK Cancel Help

26/41

FIG. 26

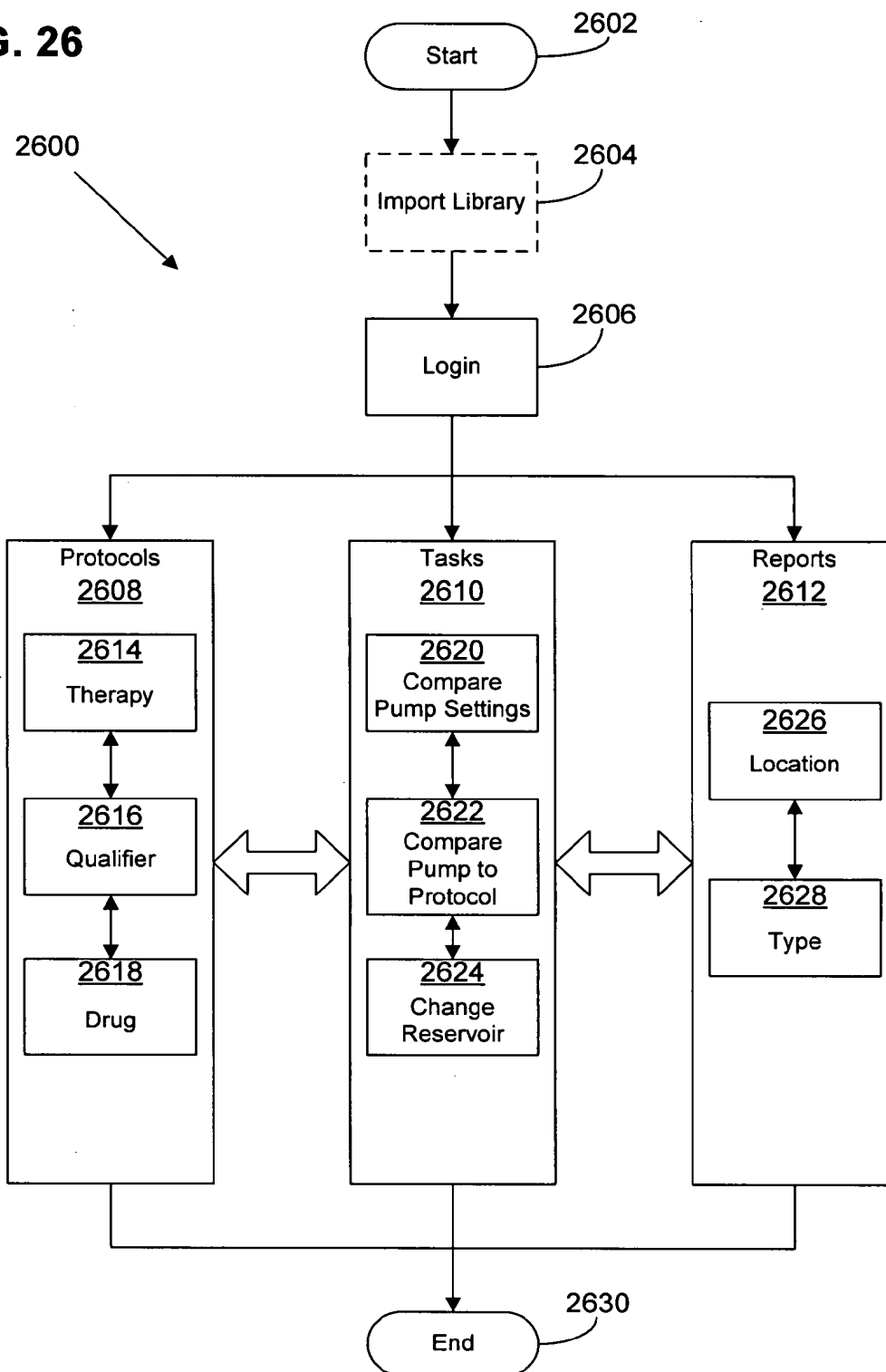


FIG. 27

Import Library to Local Computer

Help

Location

Import From Location: 2702

L:\ExportFromAdmin.mdb Browse

Library

Name: 2704

Description:

Last Modified On:

Enter Library Password:


2706 OK Cancel

2700

2708a

2708b

FIG. 28

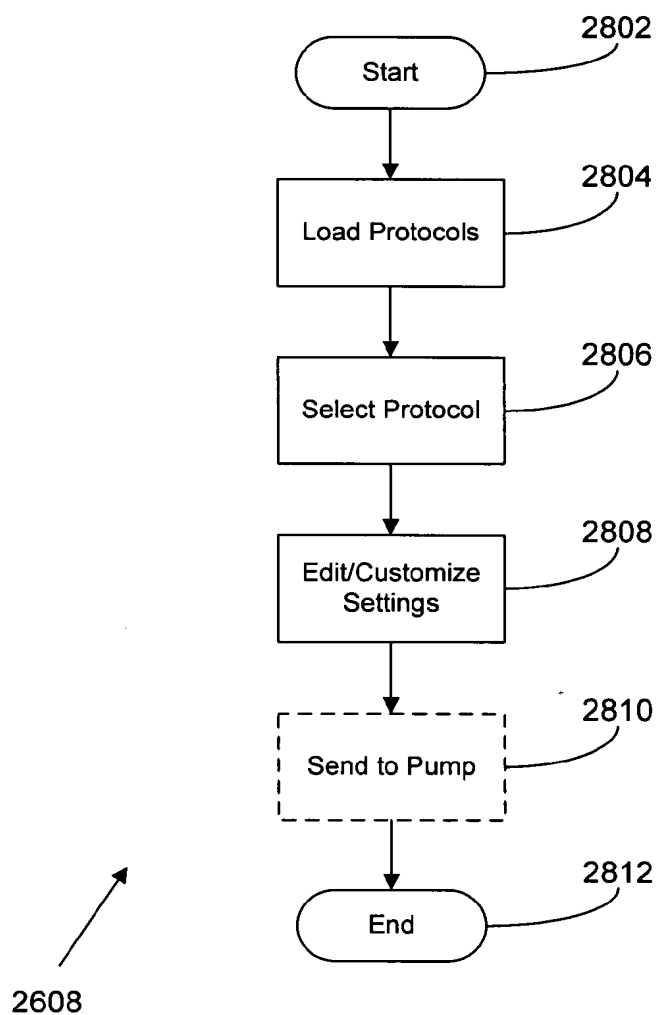

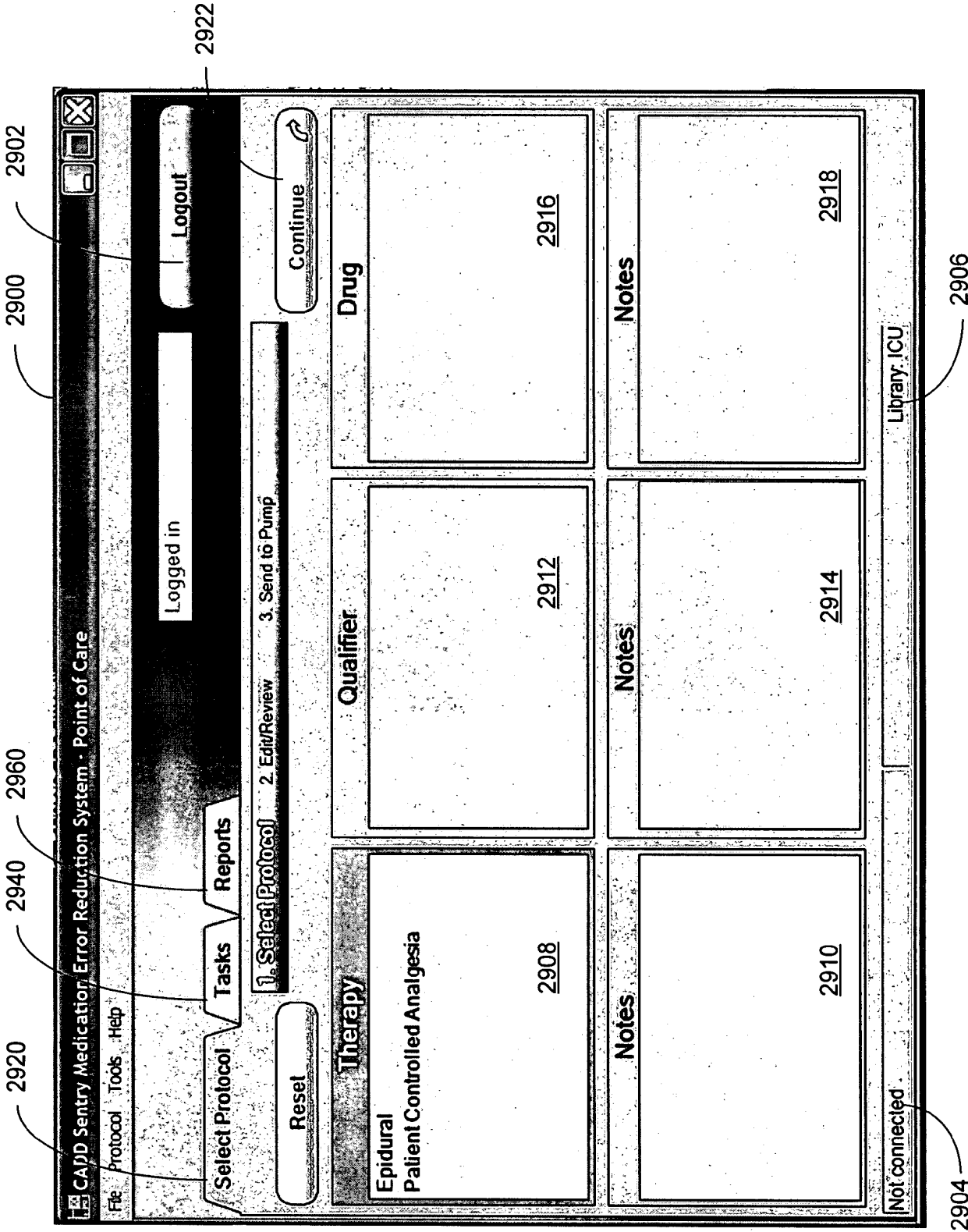
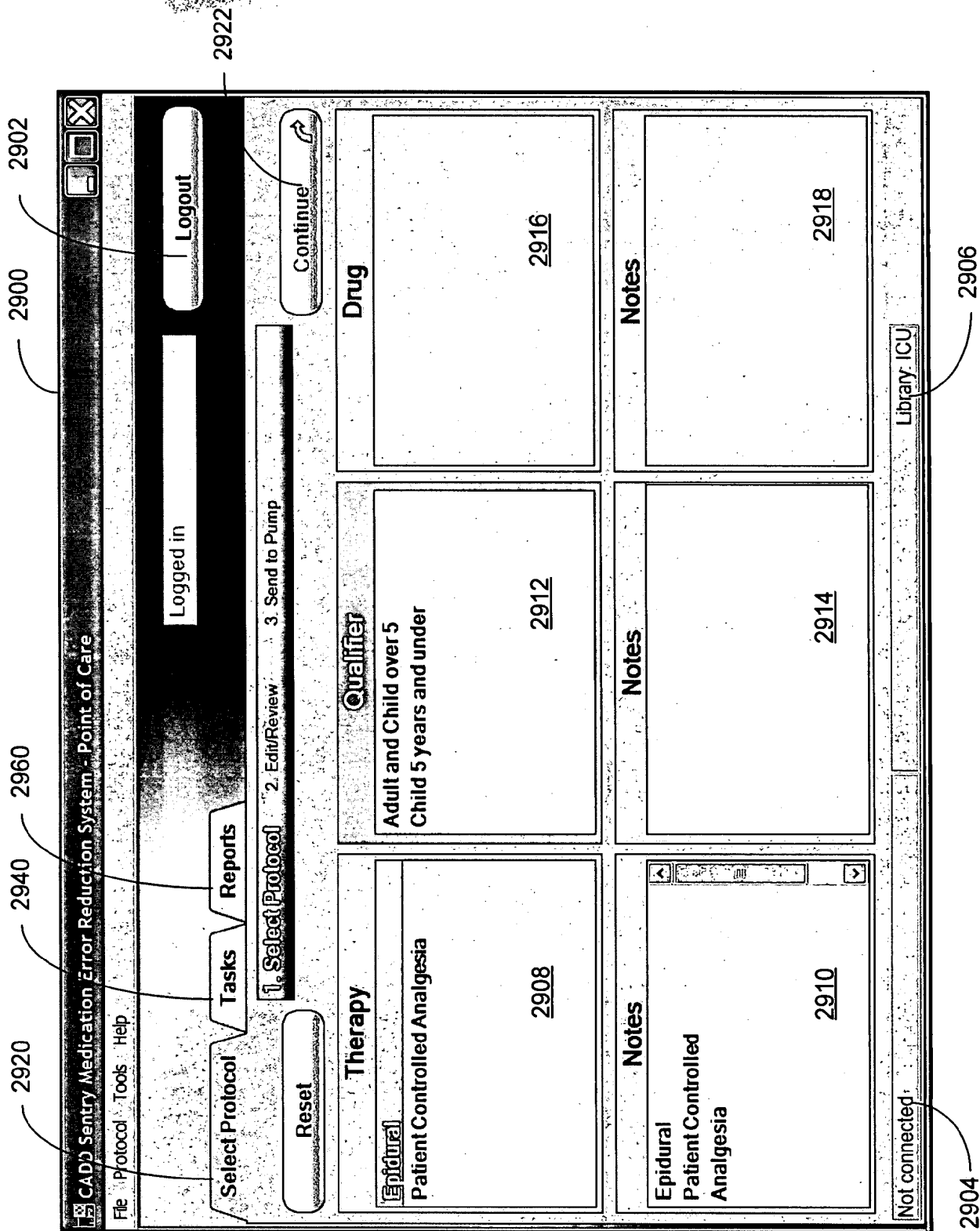


FIG. 29

FIG. 30



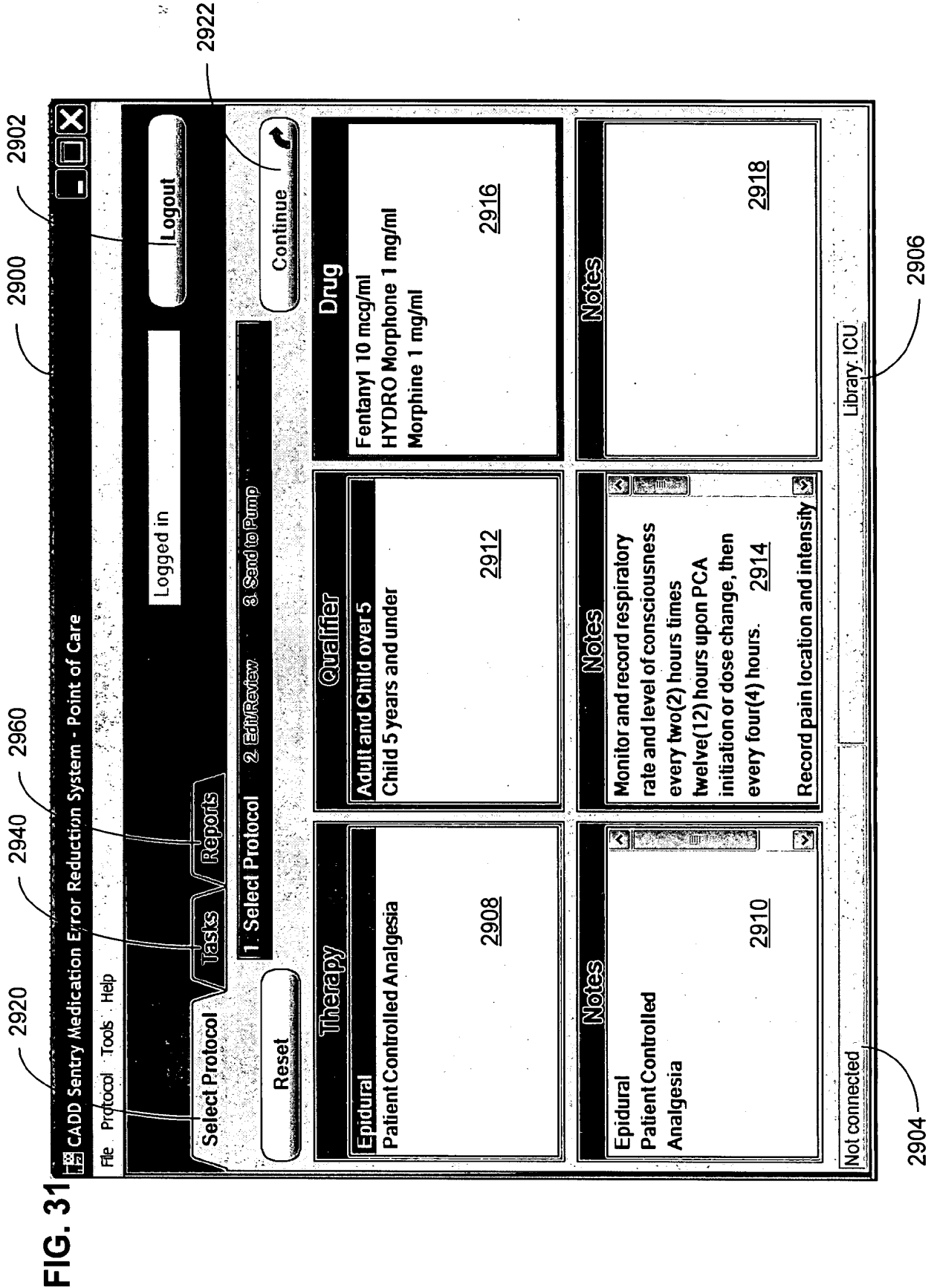


FIG. 32

CADD Sentry Medication Error Reduction System - Point of Care

File Protocol Tools Help

Logged in

Logout

1. Select Protocol 2. Edit/Review 3. Send to Pump

Reset

Therapy

Epidural

Patient Controlled Analgesia

2908

Qualifier

Adult and Child over 5

Adults over 65 years

Child 5 years and under

2912

Drug

Fentanyl 10 mcg/ml

HYDRO Morphine 1 mg/ml

Morphine 1 mg/ml

2916

Notes

Patient Controlled Analgesia

2910

Consider reducing dose gradually. Monitor for side effects such as: dizziness, SOB, anxiety, confusion, decreased urine output, hallucinations, constipation, dry mouth, nausea and vomiting. Monitor for symptoms of overdose: cold clammy

2918

Monitor and record respiratory rate and level of consciousness every two(2) hours times twelve(12) hours upon PCA initiation or dose change, then every four(4) hours. 2914

Record pain location and intensity

Not connected

Library: ICU

FIG. 33

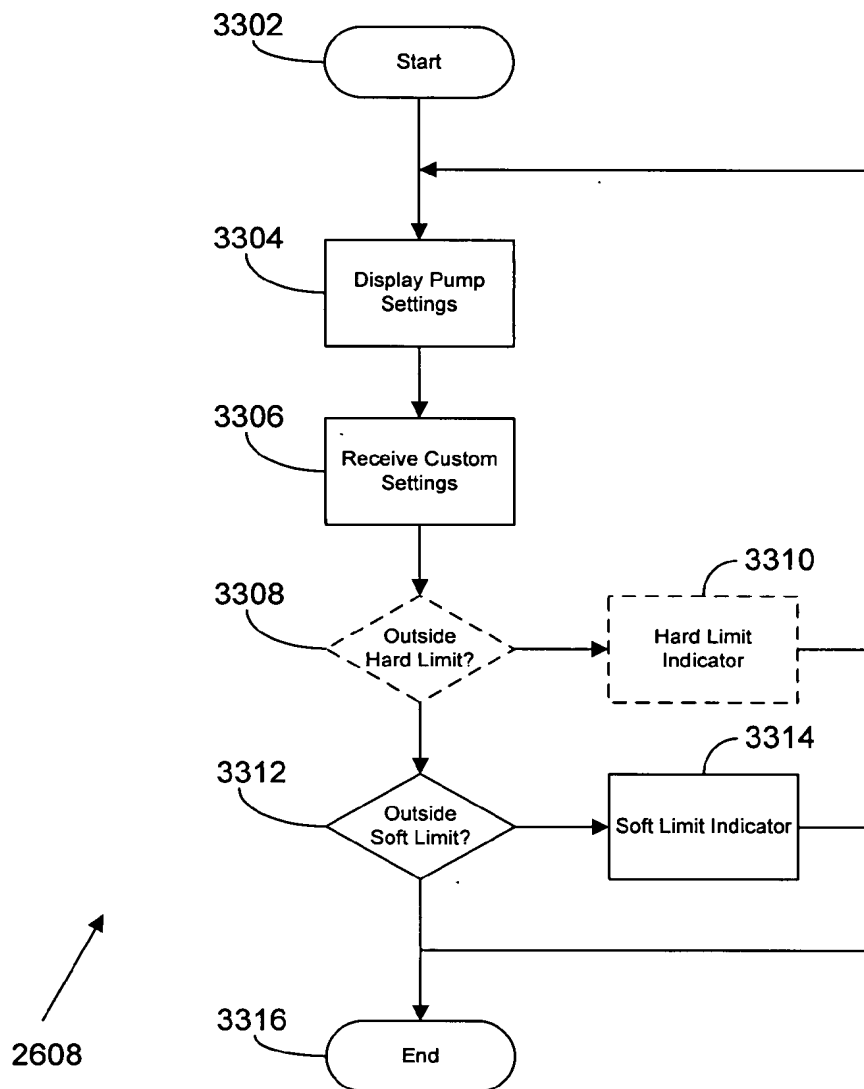


FIG. 34

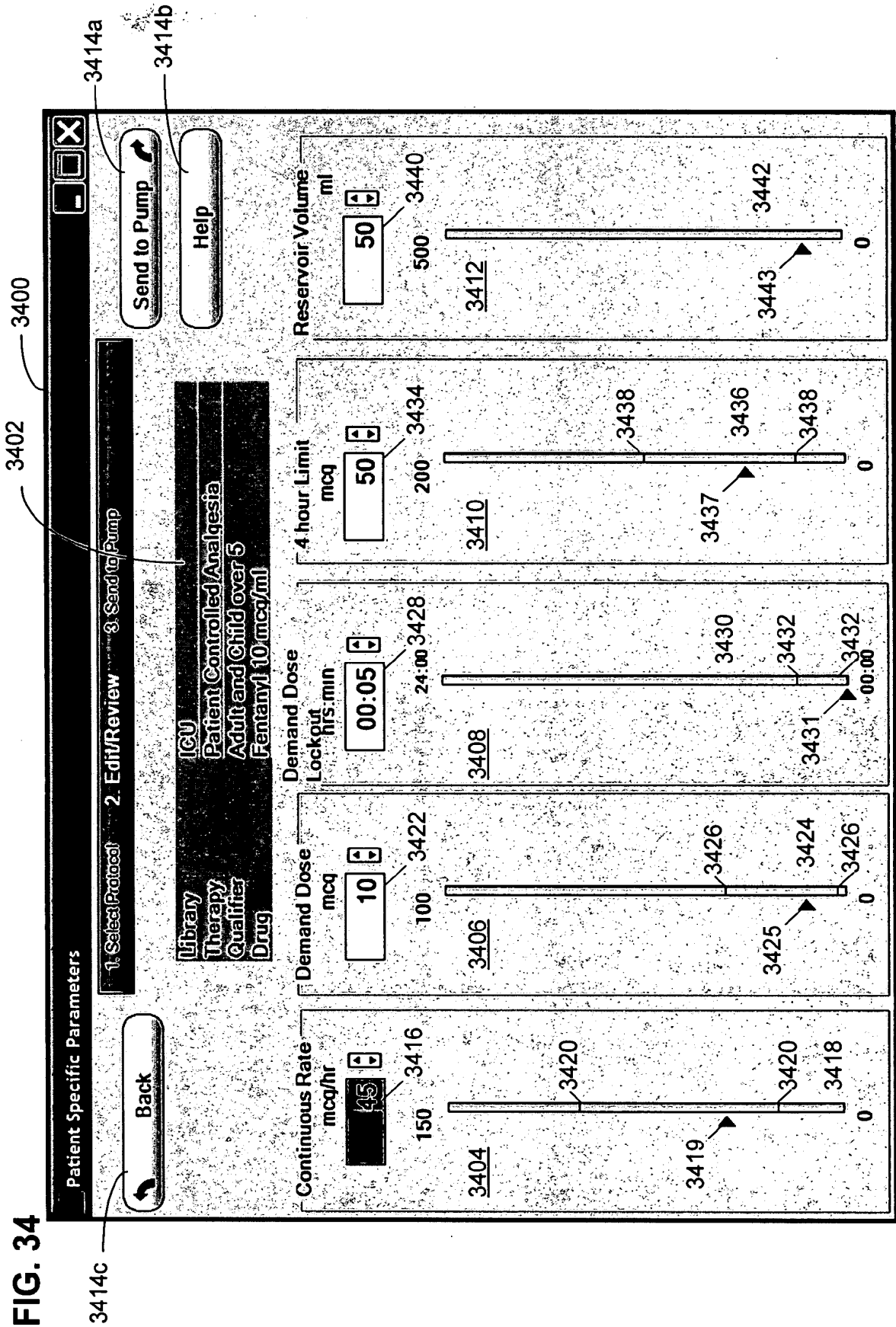


FIG. 35

Patient Specific Parameters (Modified)

1. Select Protocol 2. Edit/Review 3. Send to Pump

Back Send to Pump Help

Library ICU Epidural
Therapy Adult and Child over 5
Qualifier Fentanyl 10 mcg/ml
Drug

Continuous Rate mcg/hr
124 3416 3420 3420 3418
3404 3419 0

Demand Dose mcg
15 3422 3426 3424 3426
3406 0

Demand Dose Lockout hrs:min
00:15 3428 3430 3432 3432
3408 00:00

4 hour Limit mcg
144 3434 3438 3435 3438
3410 3437 3502 0

Reservoir Volume ml
50 3440 3442 3443 0
500

FIG. 36

Send to Pump

3604

3600

1. Select Protocol 2. Edit/Review 3. Send to Pump

Library ICU
Therapy Patient Controlled Analgesia
Qualifier Adult and Child over 5
Drug Fentanyl 10 mcg/ml

Back

Continuous Rate 45 mcg/hr 3602a

Demand 10 mcg 3602b

Demand Dose Lockout 0.5 hrs/min 3602c

4 hour Limit 50 mc 3602d

Reservoir Volume 50 ml 3602e

Please connect a pump. If a pump is connected you may need to change the COM setting.

Send 3606

Close Help

Not Connected Library ICU

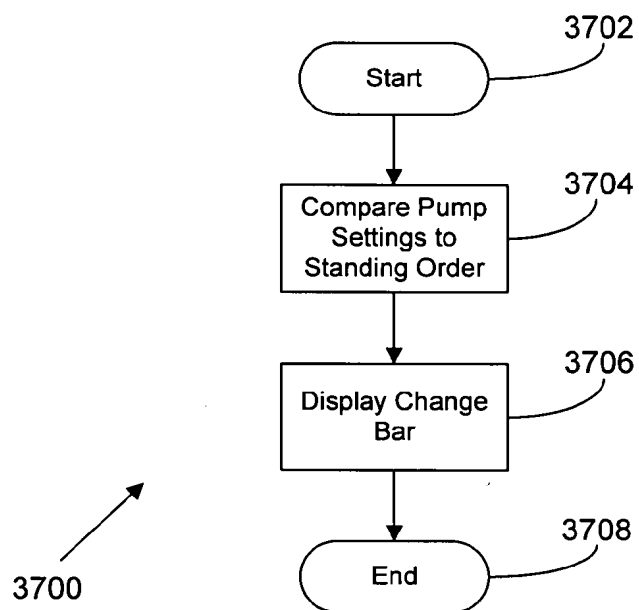
FIG. 37

FIG. 38

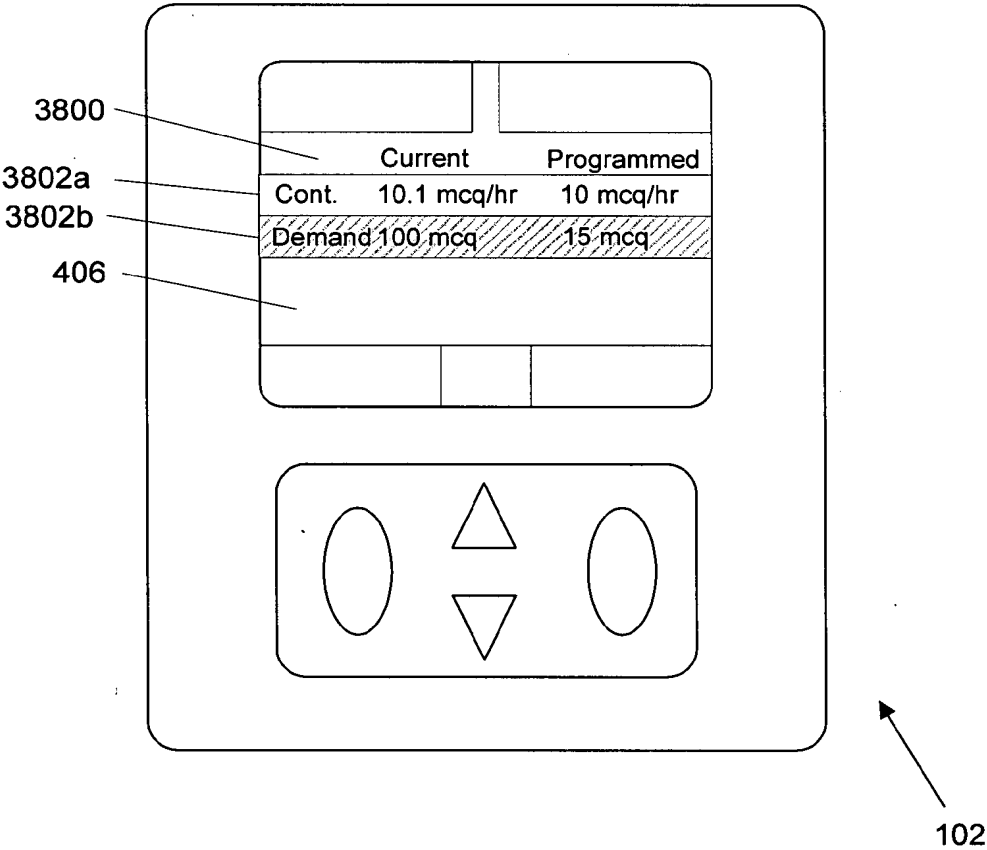


FIG. 39

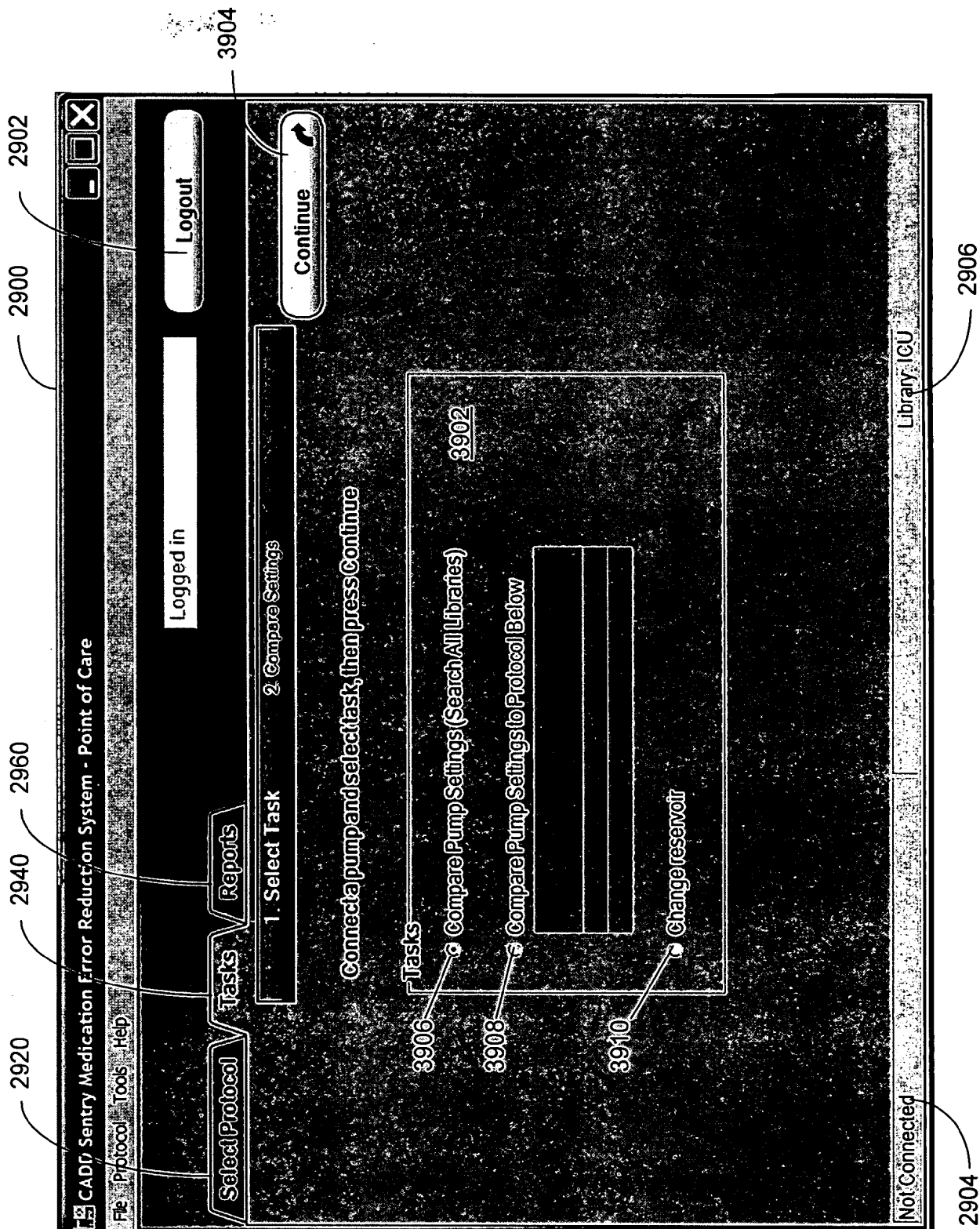


FIG. 40

